

HORIZON 6 Flex FA

Operator's Manual



 **Drucker Diagnostics**

PRODUCT FAMILY HORIZON SERIES (HORIZON 6 Flex, 6 Flex FA)

SETTINGS • CONFIGURACIÓN • PARAMÈTRES • EINSTELLUNGEN • IMPOSTAZIONI • CONFIGURAÇÕES

(EN) The LED indicator light is on for the cycle currently selected.

(ES) La luz indicadora LED está encendida para el ciclo actualmente seleccionado.

(FR) Le voyant lumineux à LED est allumé pour le programme actuellement sélectionné.

(DE) Die LED-Anzeige leuchtet für den aktuell gewählten Zyklus.

(IT) La spia LED è accesa per il ciclo attualmente selezionato.

(PT) A luz indicadora do LED está acesa para o ciclo atualmente selecionado.

	1 - Chemistry	2 - Coag (PPP)	3 - Urine
RPM (TPM U/min Giri/min)	3900 @ 10 min	3800 @ 15 Min	2200 @ 5 min
TIME (Tiempo Durée Zeit Durata Tempo)			
75-100 MM	1600 g	1500 g	500 g
125 MM	1850 g	---	600 g

SPECIFICATIONS • ESPECIFICACIONES • CARACTÉRISTIQUES • ANFORDERUNGEN • SPECIFICHE • ESPECIFICAÇÕES

	For Centrifuge Serial Number containing: YYMM04AA001	For Centrifuge Serial Number containing: YYMM41AA001
CAPACITY (Capacidad Capacité Röhrchenkapazität Capacità Capacidade)	6 (75-100 mm, 3-10 mL)	
DIMENSIONS (Dimensiones Maße Dimensioni Dimensões)	12 in x 14 in x 9 in (30 cm x 36 cm x 23 cm)	
WEIGHT (Peso Poids Gewicht Peso)	12 lbs (5.4 kg)	
SOUND LEVEL (Nivel de ruido Niveau sonore Lautstärkepegel Rumorosità Nivel de ruido)	64 dB A	
ROTOR RADIUS (Radio Rayon Raggio Raio)	3.75 in (9.5 cm), 4.25 in (11 cm)	
CENTRIFUGE MOTOR (Moteur Motore)	½ H.P. Brushless DC	
CYCLE TIME (Tiempo del ciclo Durée du cycle Zykluslaufzeit Durata ciclo Tempo de ciclo)	1 – 30 min (+/- 2%)	
SUPPLY VOLTAGE (Voltaje Tension Spannung Tensione Tensão)	100 – 240 (+/- 10%)	
SUPPLY FREQUENCY (Frecuencia Fréquence Frequenz Frequenza Frequência)	50/60 Hz	
CURRENT CONSUMPTION (Energía necesaria Alimentation Energiebedarf Requisiti di alimentazione Requisito de potência)	220 Watts	2.2A at 115 VAC 1.1A at 230 VAC
ENVIRONMENTAL RANGE (Rango de temperatura ambiental Plage de température de fonctionnement Umgebungsbedingungen Intervall o temp. ambiente Faixa ambiental)	16 – 32° C	
MAX G-FORCE (Fuerza g máxima Force G maximale Maximale g-Kraft Forza G)	1600 g / 1850 g	

massima Força G máxima)	
MAX SPEED (Velocidad máxima Vitesse maximale Maximale Beschleunigung Velocità massima Velocidade máxima)	3900 RPM

ENVIRONMENTAL CONDITIONS • CONDICIONES MEDIOAMBIENTALES • CONDITIONS ENVIRONNEMENTALES • UMWELTBEDINGUNGEN • CONDIZIONI AMBIENTALI • CONDIÇÕES AMBIENTAIS

SET-UP SITE (Sitio de instalación Site de mise en place Standorteinrichtung Sito di installazione Local de montagem)	(EN) Indoor Use Only (ES) Solo para uso en interiores (FR) Uniquement pour un usage à l'intérieur (DE) Nur für den Innenbereich (IT) Solo per uso interno (PT) Apenas para uso interno
ALTITUDE (Altitud Höhenlage Altitudine)	(EN) Up to 2,000m from Sea Level (ES) Hasta 2000 m sobre el nivel del mar (FR) Jusqu'à 2 000 m au-dessus du niveau de la mer (DE) Bis zu 2.000 m über dem Meeresspiegel (IT) Fino a 2.000 m sul livello del mare (PT) Até 2.000 m do nível do mar
AMBIENT TEMPERATURE (Temperatura ambiente Température ambiante Umgebungstemperatur Temperatura do ambiente)	(EN) 5 °C to 40 °C (ES) De 5 °C a 40 °C (FR) 5 °C à 40 °C (DE) 5 °C bis 40 °C (IT) Da 5 °C a 40 °C (PT) 5 °C a 40 °C
HUMIDITY (Humedad Humidité Luftfeuchtigkeit Umidità Umidade)	(EN) Maximum relative humidity 80% for temperatures up to 31 °C, decreasing linearly to 50% relative humidity at 40 °C (ES) Humedad relativa máxima del 80 % para temperaturas de hasta 31 °C, disminuyendo linealmente hasta el 50 % de humedad relativa a 40 °C (FR) Humidité relative maximale de 80 % pour des températures pouvant atteindre jusqu'à 31 °C, Décroissance linéaire jusqu'à 50 % d'humidité relative à 40 °C (DE) Maximale relative Luftfeuchtigkeit 80 % bei Temperaturen bis zu 31 °C, linear abnehmend auf 50 % relative Luftfeuchtigkeit bei 40 °C (IT) 80% di umidità relativa massima per temperature fino a 31 °C, decrescente in modo lineare fino al 50% di umidità relativa a 40 °C (PT) Umidade relativa máxima de 80% para temperaturas de até 31 °C, diminuindo linearmente para 50% de umidade relativa a 40 °C
OVERVOLTAGE CATEGORY (Categoría de sobretensión Catégorie de surtension Überspannungskategorie Categoria di sovratensione Categoria de sobretensão)	II
POLLUTION DEGREE (Grado de contaminación Degré de pollution Verschmutzungsgrad Grado di inquinamento Grau de poluição)	2

Contents

- English 1
- Symbols..... 2
- Caution and Warning Statements 3
- Model Description 5
 - Features 5
 - Intended Use..... 5
 - Warranty 5
- Initial Setup..... 6
- Quick Start 7
 - Single Cycle Lock 7
 - Preset Lock 8
- Setting or Modifying a Cycle 8
- Display Cycle Count 9
- Operation..... 9
- Balancing Loads 9
- Care and Preventative Maintenance 10
- Cleaning and Disinfection 10
- Troubleshooting..... 11
- Calculating the G-Force 13









6 FLEX FA

English

Operator's Manual



SYMBOLS

Symbol	Definition	Use
	Caution	Caution to safety hazard. Potential risk of personal injury or damage to the instrument if improperly handled. Consult the manual before proceeding.
	Manufacturer	Manufacturer of record.
	Electrical and electronic products recycling symbol	Recycle only as electronic waste. Do not dispose in normal waste.
	RoHS Compliant	Compliance with RoHS environmental standards.
	CE Mark	Denotes conformity to specific European directives and regulations.
	UK Mark	Denotes conformity to specific UK directives and regulations.
	MET Listing	Denotes conformity to specific safety standards and regulations.
FDA LISTED	FDA Listed	Denotes that the product has been properly listed with the FDA.
	ISO Certification	Denotes conformity to quality standards and quality management systems.

CAUTION AND WARNING STATEMENTS



This device is intended to be operated by properly trained personnel who have carefully read the operating manual and are familiar with the function of the device. Users should also comply with the specimen receptacle manufacturer's specific instructions for use, in addition to any other protocols established by the testing organization.



WARNING: For the safety of both the operator and service personnel, care should be taken when using this centrifuge if handling substances that are known to be toxic, radioactive or contaminated with pathogenic microorganisms. Use appropriate personal protection equipment (PPE). When Risk Group II materials are used, (as identified in the World Health Organization "Laboratory Bio-Safety Manual"), a Bio-Seal should be employed. In the event that materials of a higher risk group are being used, more than one level of protection must be provided. The use of flammable or explosive materials as well as those materials which have a vigorous chemical reaction is prohibited.



Unplug the centrifuge before cleaning or performing maintenance.



WARNING: Inspect centrifuge for cracks or physical damage to cabinet, lid, rotor, or tube holders. Damage may result in unsafe operation. Discontinue use until repairs have been performed.



The use of flammable or explosive materials as well as those materials which have a vigorous chemical reaction is prohibited.



For your safety and durability of the machine, never transport or store centrifuge with tube holders inside the machine.



WARNING: "Universal precautions"¹ should be followed in handling all items contaminated with blood or other bodily fluids.



This equipment generates, uses, and can radiate radio frequency energy, and if not installed and used in accordance with this operator manual, may cause interference to radio communications.



Operation of this equipment in a residential area may cause interference, in which case the user will be required to correct the interference at his own expense.



Operation of this equipment in a manner not specified by the manufacturer may impair the protection provided by the equipment.



Electrical Safety protection is provided by properly connecting the centrifuge to earth ground. Use only the manufacturer provided line cord and ensure that it is connected to a properly grounded power receptacle. Failure to do so will result in an electrical hazard.



WARNING: Do not make modifications to or remove any hardware from rotor without prior authorization from Drucker Diagnostics.



WARNING: Only use Drucker Diagnostics components in this centrifuge.



Due to the lack of the possibility of human exposure, all Drucker centrifuges and accessories sold by Drucker Diagnostics, Inc. are compliant without any special labeling required by the California Safe Drinking Water and Toxic Enforcement Act (Proposition 65).

1 Recommendations for Prevention of HIV Transmission in Health Care Settings. MMWR 1987; 36 (Supplement #2S)

MODEL DESCRIPTION

HORIZON is a versatile line of centrifuges pre-programmed with the most convenient cycle settings for biological sample processing. Cycle settings can be changed to accommodate custom settings.

This general purpose laboratory centrifuge may also be used to spin approved containers with biologics, chemicals (non-flammable, non-explosive, non-volatile, and non-highly reactive), and environmental samples.

FEATURES

- The first three (3) cycles are conveniently pre-set and labeled for your most common applications. Use the default cycles or customize them as needed.
- If desired, the control panel can be temporarily locked on one cycle for error-free reproducibility.
- Up to 10 cycles can be programmed for time, speed, and braking and labeled with a custom name. Cycles can be programmed by g-force (RCF) rather than speed to facilitate matching validated cycles and manufacturers IFUs.
- Lid lighting indicates the centrifuge's status (ready, running, done), informing the operator when tubes are ready for the analyzer and preventing tubes from being left in the centrifuge longer than necessary (patent pending).
- A traditional audible alert indicates the completion of the cycle. The audible alert can be muted.
- Cool-Flow design prevents overheating of samples by using ambient air to keep specimens at room temperature.
- The tube holders are fiber reinforced for high strength, durability, and years of trouble-free use.
- A clear lid permits safe observation of samples and optical calibration of speed.
- The lid safety system prevents the centrifuge from operating unless the lid is closed and latched.
- The lid safety system only allows entry into the centrifuge after the rotor has completely stopped.
- The high-power brushless motor provides years of operation with no routine maintenance.

INTENDED USE

General purpose laboratory centrifuge, intended for the density-based separation of fluids through centripetal acceleration.

WARRANTY

Drucker Diagnostics warranties that this centrifuge is free from defects in workmanship and parts for 2 years.


INITIAL SETUP

- Unpack and verify that all the following are included:
 - Centrifuge
 - Power cord
 - Tube holders
 - Quick Start Insert
- Setup the centrifuge on flat and level surface. A bench top clearance height of 21" (54 cm) is required to open the lid.
- The centrifuge should have 6" (15 cm) of clear space around the centrifuge. Proper ventilation is necessary to prevent the overheating of samples as well as premature failure of the centrifuge. Choose an area which allows unencumbered air flow, and where the temperature remains between 16°C and 32°C.
- No hazardous material shall be permitted in the clearance envelope during operation.
- The operator time within the envelope shall be limited to the time necessary for loading, unloading, and centrifuge operation only.
- Plug the line cord into the centrifuge.
- Plug the line cord into an approved electrical outlet.
- Turn on the power switch in the back of the centrifuge.



BE SURE THE ELECTRICAL OUTLET IS ALWAYS ACCESSIBLE AS THE LINE CORD IS THE MEANS OF EMERGENCY DISCONNECTION!

QUICK START

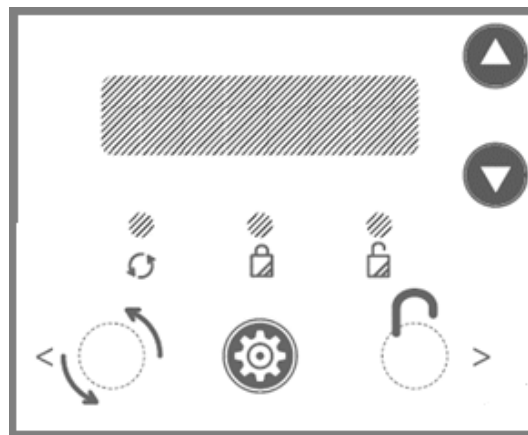




	Start	Begins running the cycle displayed on the screen. The lid must be closed.
	Unlock	Allows access into the rotor chamber by engaging the unlocking mechanism. Entry is only possible when the rotor is stopped.
	Stop	Pressing the UNLOCK button during operation will terminate the run and unlock the lid after the rotor has come to a stop.
	Cycle Selection	Press the up and down buttons next to the screen to select the desired saved cycle.

To ensure repeatability, the centrifuge can be locked either on one cycle (Single Cycle Lock) or restricted to the saved cycles (Preset Lock). The Single Cycle Lock also prevents making changes to the parameters of the selected cycle. The Preset Lock allows selection of any saved cycle, but prevents changing the parameters of saved cycles.

SINGLE CYCLE LOCK








	Lock on Single Cycle	Select desired cycle. Press and hold the UNLOCK button for 5 seconds. Two beeps will confirm that cycle selection is locked.
	Cancel Single Cycle Lock	To re-enable cycle selection, press and hold the UNLOCK button for 5 seconds. Three beeps will confirm that cycle selection is now unlocked.

PRESET LOCK

	Enter Settings Mode	Press the GEAR button.
	Lock Saved Cycles	Press and hold the UNLOCK button for 5 seconds. Two beeps will confirm that the Preset Lock is active.
	Cancel Preset Lock	With lid open and while not in Settings mode, press and hold the UNLOCK button for 5 seconds until two beeps are heard. Saved cycles can again be edited or programmed.




SETTING OR MODIFYING A CYCLE

Up to 10 cycles can be named, programmed for time, speed, and braking, and saved in memory.

	Cycle Selection	Press the up and down buttons next to the screen to select the desired cycle.
	Enter Programming Mode	Press the GEAR button.
	Navigate between Parameters	Use the left and right arrow button until the parameter to be edited is blinking and underlined.
	Change Parameter Values	To change the parameter shown on the display, use the up and down buttons next to the screen.
	Displaying the Correct g-force	HORIZON 6 Flex FA only: To display the real g-force (RCF) for the cycle, the right tube holder must be selected. Universal blue tube holders accommodate 75 and 100 mm tubes. 125 mm tubes require the orange tube holder. If the tube holder selection in the program does not match the tube holder used, the RCF number displayed will not reflect the actual value.
	Naming the Cycle	Navigate to the cycle name with the left and right arrow. Change the blinking character of the name with the up and down buttons, then move to the next character with the right arrow. Repeat.
	Save and Exit Programming Mode	Press the GEAR button. The cycle displayed on the screen will be automatically saved.

DISPLAY CYCLE COUNT

These actions will display the cycle count.

	Enter Settings Mode	Press the GEAR button.
	Navigate to Cycle Counter	Press the left arrow once to access the cycle counter
	Exit Settings Mode	Press the GEAR button.

OPERATION

- Place the tubes into the tube holders. Be sure to follow the rules for balanced loads as listed in the next section.
- Close the lid and turn the lid knob clockwise to its complete stop position.
- The digital screen shows the currently selected cycle. To select another cycle, press the UP or DOWN button in succession until the desired cycle is selected.
- Pushing the START button on the control panel will start the spin cycle.
- When the cycle is completed, the rotor will slow to a complete stop and the lid light will flash.
- The unlocking mechanism will engage for 60 seconds allowing entry into the rotor chamber. To unlock after more than 60 seconds have elapsed, press the UNLOCK button. The lid will unlock for another 15 seconds.
- Turn the lid knob counterclockwise and open the lid. The lid light will turn off.


BALANCING LOADS



Your centrifuge must contain a balanced load to work properly. Spinning balanced loads will extend the life of the centrifuge and produce better results. Use the following rules when loading the rotor. If an odd number of samples is to be spun, fill a tube with water to match the weight of the unpaired sample and place it across from this sample.

*Opposing tube holders must be equally loaded or empty or loaded with equally weighted samples.
When loading only 3 tubes, they must be of equal weight.*

6 Tube Centrifuges

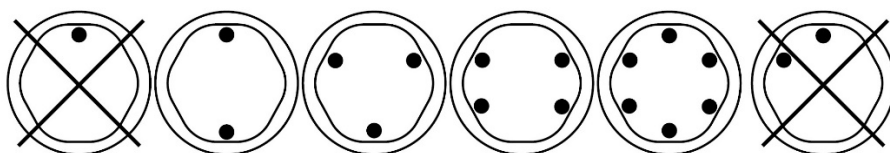

CARE AND PREVENTATIVE MAINTENANCE

With proper care and maintenance, your centrifuge will provide years of laboratory service. For proper care, the following steps should be taken:

- **Always Spin Balanced Loads:** Make certain that you are always spinning a balanced load, as shown in the previous section. These centrifuges have a unique counter balanced motor mounting design which produces excellent vibration dampening. However, out-of-balance loads may break glass test tubes and may produce unsatisfactory separation results. Proper load balancing will improve sample separation and extend the life of the centrifuge.
- **Motor and Electrical Maintenance:** The highest quality electrical components have been selected for the centrifuges and should not need maintenance or servicing for the life of the centrifuge.
- **Tube Holder Replacement:** It is recommended that the tube holders be replaced after 24 months of use. Inspect tube holders regularly for cracks. If cracks are discovered, replace immediately.
- **Remove Accessories Before Moving:** All tube holders, samples, and caps must be removed from the rotor chamber before transporting or storing the centrifuge to prevent damage and injury.

CLEANING AND DISINFECTION

To prolong the life of the centrifuge, cleaning and disinfection is recommended every six months or whenever there is a spillage or tube breakage. Contaminants must be removed immediately or corrosion and premature degradation of components can occur. Before using any cleaning or decontamination methods other than those recommended by the manufacturer, users should verify with the manufacturer that the proposed method will not damage the equipment.



Cleaning and Decontamination may be necessary as a safeguard before laboratory centrifuges, rotors, and any accessories are maintained, repaired, or transferred.

- Unplug the centrifuge before cleaning.
- Use appropriate personal protective equipment (PPE).
- Apply Cleaning solutions with a dampened towel or cloth ONLY. Do not spray or pour cleaning solution directly onto or into the centrifuge. Do not saturate or submerge the centrifuge in water or other cleaning solutions as this will cause damage, create a safety risk, and void the warranty.
- ONLY isopropyl alcohol or a 10% (5500 PPM) bleach solution should be used to disinfect the centrifuge and its accessories.



All surfaces must be dried immediately after cleaning and disinfecting.


TBQ GERMICIDAL PRODUCTS ARE NOT RECOMMENDED AS THEY MAY CAUSE DAMAGE TO THE CENTRIFUGE. REFRAIN FROM USING TO PREVENT VOIDING THE WARRANTY.

- Fully/partially halogenated hydrocarbons, ketones, esters, ethers, benzyls, ethyl benzenes, and all other chemicals not prescribed by the manufacturer shall not be used as they may cause damage to the rotor chamber, rotor, tube holders, accessories and centrifuge exterior and void the warranty.

TROUBLESHOOTING

NOTE: The latch must be turned completely clockwise to its stop position for the centrifuge to operate.

<p>The centrifuge does not run</p>	<ul style="list-style-type: none"> ○ Verify that the centrifuge is powered. The screen should be lit. ○ If “Lid not closed” message is displayed, make sure the lid latch is turned completely clockwise to its stop position. ○ If the centrifuge still does not run, contact Customer Service.
<p>The rotor does not spin freely</p>	<ul style="list-style-type: none"> ○ Make sure nothing has fallen into the rotor chamber, following the procedure above. ○ If nothing obstructs the rotor, the rotor may be damaged. Contact Customer Service for further assistance.
<p>The centrifuge makes a rattling noise when running</p>	<ul style="list-style-type: none"> ○ Stop the centrifuge. Open the lid. ○ Wearing PPE, remove tubes and tube holders/buckets and look for fallen objects or debris. Carefully reach inside the rotor chamber with a tool to remove them. ○ Inspect the rotor, tube holders or buckets for damage. ○ If the tube holders or buckets have any damage, even slight, safely dispose of them and replace them. ○ If the rotor appears damaged, contact Customer Service for further assistance.
<p>Excessive noise or vibration when the centrifuge is running</p>	<ul style="list-style-type: none"> ○ Verify that all four centrifuge feet are properly seated on a flat surface. ○ Ensure that the load is balanced according to instructions in the “Balancing Loads” section of this manual. ○ Make sure that nothing has fallen into the rotor chamber.
<p>The centrifuge stops and beeps continuously</p>	<p>The load is not balanced. Press the UNLOCK button, open the lid, and balance the load as recommended elsewhere in this manual.</p>
<p>The centrifuge is stuck on one of the settings</p>	<p>Cycle selection is locked. Press the UNLOCK button for 5 seconds.</p>
<p>The cycle time and speed are not set to the desired value</p>	<p>Check the setting by following the instructions in the section on Changing Cycle Settings. If the preset is not the desired length, follow the procedure on the same page to change the preset time.</p>
<p>Cycle parameters cannot be changed</p>	<ul style="list-style-type: none"> ○ If cycle selection is locked on one cycle, press the UNLOCK button for 5 seconds. Then, press the GEAR button and follow the instructions elsewhere in this manual. ○ If different saved cycles can be selected but not modified, the centrifuge is in Preset Lock mode. Press the GEAR button, then the UNLOCK button for 5 seconds. You should now be able to change cycle parameters.
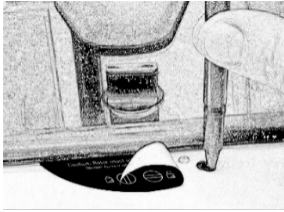
<p>The centrifuge does not unlock after a run is completed</p>	<ul style="list-style-type: none"> ○ Wait until the rotor has come to a complete stop. If the lid knob still cannot be rotated, press the UNLOCK button and try again. <ul style="list-style-type: none"> ○ If no LED light is on, the unit is not powered and the lid will not unlock by conventional means. Remove the latch label and use a pen to manually disengage the locking mechanism. Pull the mechanism towards the control panel and then unlatch and open the lid. ○ If the unit is damaged, contact Customer Service for assistance. 
<p>The lid does not open</p>	<ul style="list-style-type: none"> ○ Ensure that the lid knob is turned fully counterclockwise. ○ If the knob cannot be turned counterclockwise, turn it fully clockwise, press UNLOCK, and turn counterclockwise. ○ If the lid remains locked after this and will not unlock, the electronics may have been damaged. Contact customer service for assistance.
<p>Clicking noise during braking gets loud</p>	<ul style="list-style-type: none"> ○ Make sure that the screw in the center of the rotor is tight.
<p>Lid does not stay up</p>	<ul style="list-style-type: none"> ○ Tighten the center screw on the lid hinge.

CALCULATING THE G-FORCE

The I.F.U.s of tube manufacturers recommend cycles at a minimum G-Force, which can be calculated if you know the RPM and the radius. Use the formula below.

In Centimeters:
$$\text{RCF or G-force} = 0.00001118 \times \text{Rotor Radius (cm)} \times (\text{RPM})^2$$

In Inches:
$$\text{RCF or G-force} = 0.0000284 \times \text{Rotor Radius (in)} \times (\text{RPM})^2$$

Complies with UL61010-1/CSA C22.2 No. 61010-1 and IEC61010-2-020

Protected by U.S. Patents #6,811,531, # 7,422,554, #D718,463, & #D734,489. Other Patents Pending

INSTRUCTIONS FOR DISPOSAL OF WEEE BY USERS IN THE EUROPEAN UNION



This product must not be disposed of with other waste. Instead, it is the user's responsibility to dispose of their waste equipment by handing it over to a designated collection point for the recycling of waste electrical and electronic equipment. The separate collection and recycling of your waste equipment at the time of disposal will help to conserve natural resources and ensure that it is recycled in a manner that protects human health and the environment. For more information about where you can drop off your waste equipment for recycling, please contact your local city office, waste disposal service, or where you purchased the product.

PRODUCT FAMILY HORIZON Series (HORIZON 6, 6FA, 12, 24)

MANUAL 03-0-0002-0169, Rev. D

**FDA
LISTED**

