



PORTABLE LED SURGICAL HEADLIGHTS

Superior Surgical Lighting

The Color Rendition Index (CRI) is a quantitative measure of a light's ability to reveal full spectrum colors and textures. Higher ratings indicate greater fidelity to the colors of tissues and skin and can reduce the time to diagnosis and treatment. Headlights with a high CRI are desirable in color-critical applications including OB/Gyn, Dermatology, Oral Maxillofacial Surgery, ENT, Neuro/Spine, and more.

Important Notes

- Headbands and headlight modules cannot be immersed in disinfecting liquids or sterilizers.
- Keep ends of cable clean by wiping with a damp cloth. Allow to dry before use.
- Clean the optics only with lens tissue available at any camera store. Follow the instructions on the package.
- To ensure long life for the headgear, cable and power source, store them in a cool dry environment.

**MADE
in the
USA**



Magnum

SPECIFICATIONS


Illumination	30, 45, 75k Variable Lux at 12"
Color Rendition Index (CRI)	70
Color Temperature	5,000 °K
Spot Size	10mm to 110mm
Operating Time	3.2 hrs at 100%
Recharge Time	3 hrs
Weight of Headlight	300 grams
Weight of Power Pack	350 grams
Headband	Dual Ratchet or Soft Velcro



IsoLED Plus+

SPECIFICATIONS


Illumination	45,000 Lux at 12"
Color Rendition Index (CRI)	92+
Color Temperature	4,000 °K
Spot Size	10mm to 110mm
Operating Time	4.25 hrs at 100%
Recharge Time	2.5 hrs
Weight of Headlight	290 grams
Weight of Power Pack	190 grams
Headband	Dual Ratchet or Soft Velcro

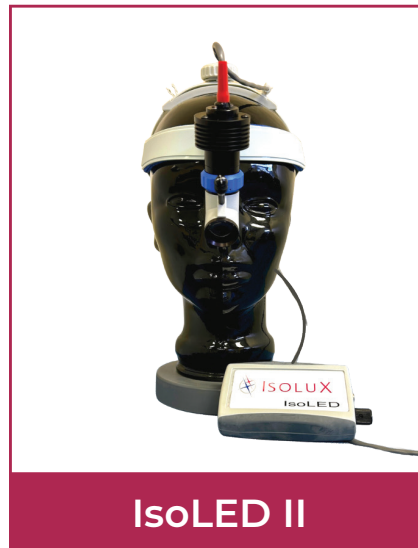

IsoLED II

SPECIFICATIONS

Illumination	25,000 Lux at 12"
Color Rendition Index (CRI)	95
Color Temperature	4,500 °K
Spot Size	10mm to 110mm
Operating Time	6 hrs at 100%
Recharge Time	3.75 hrs
Weight of Headlight	290 grams
Weight of Power Pack	190 grams
Headband	Dual Ratchet or Soft Velcro

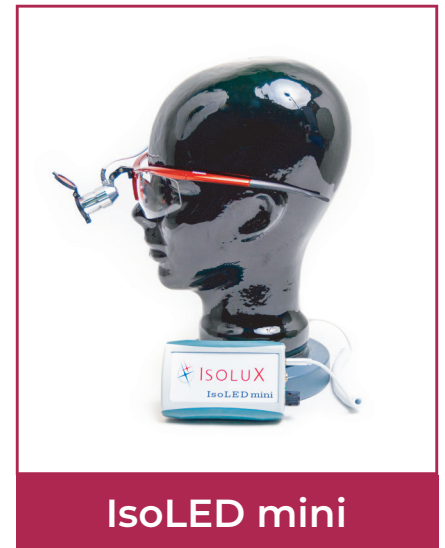

IsoLED mini

SPECIFICATIONS

Illumination	up to 22,000 Lux at 12"
Color Rendition Index (CRI)	92
Color Temperature	4,500 °K
Spot Size	60 mm @ 12"
Operating Time	10hrs typical at 100%
Recharge Time	3.6 hrs
Weight of Headlight	x3 Variations
Weight of Power Pack	190 grams
Headband Variations	Available