

THE 502KL

PRIMARY CARE SCALE WITH INTEGRATED DIGITAL HEIGHT ROD

The newest member of the number-one-selling family of physician scales, the all-digital Health o meter® 502KL features a fully integrated digital height rod and EMR capabilities, providing the optimum in quality, reliability, and ease of use. The 502KL, created to cater to various healthcare settings like doctor's offices, hospitals, clinics, and acute and long-term care facilities, showcases a market-leading built-in digital height rod, a spacious non-skid platform, and a reinforced base that ensures enhanced stability. The scale offers providers and patients an eye-level weight indicator, accommodating patients weighing up to 660 pounds. Its user-friendly digital interface provides a range of functions, including BMI calculation, hold/release, zeroing, taring, reweighing, recalling previous measurements, switching between KG and LB units, and clear and enter functions. The 502KL offers the best in its class extra-wide 8" wheelbase making transport easy. The exclusive Everlock® function allows the measuring unit (LB or KG) to be permanently locked on the scale.

The industry-leading physician scale provides easy integration and seamless communication with today's EMR applications. The 502KL offers the option to incorporate wireless connectivity for interfacing with computers, EMR software, or other electronic devices.

With its low platform height, the 502KL is easy to set up and offers patients a more comfortable and stable weighing experience.

500HB

"live" wrap-around handlebar

The 500HB is a unique handlebar accessory for the 502KL and is easy to attach to any scale using the pre-drilled mounting holes. It provides the stability that many patients need without affecting the accuracy of the weight reading.



THERE'S STRENGTH IN OUR NUMBERS.

When you add up the numbers on our scale, it's not surprising that the 500KL family of scales is the most widely used scale in physician offices nationwide.

660 LB
CAPACITY

150,000+
INSTALLED
NATIONWIDE

120+ LOCAL
FIELD SPECIALISTS

8" WIDE WHEELBASE
BEST IN ITS CLASS


#1 TRUSTED
PHYSICIAN SCALE

100 YEARS
OF SERVICE

Health o meter®
Professional

Health o meter® Professional 502KL

weigh**easier**®

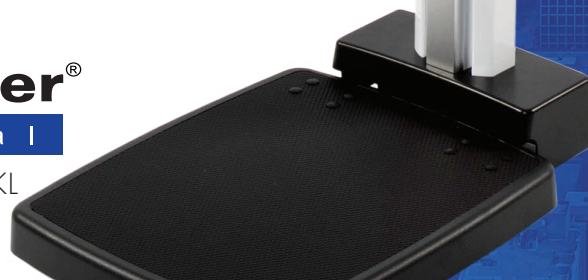

TECHNICAL SPECIFICATIONS

Capacity

660 lb / 300 kg

Resolution

0.2 lb / 0.1 kg

EMR connectivity

Via optional Pelstar® wireless technology

Display

1" / 25 mm, LCD display

Batteries

6 AA batteries included, optional adapter (order# ADPT31)

Functions

LB/KG conversion
LB/KG lock out
BMI
Zero
Tare, Pre-tare
Hold/Release
Auto Zero
Auto Off
Reweigh
Recall
Parent/Child weighing function
Everlock® - permanent unit lock

Platform Size (w x d x h)

13 3/4" x 16 1/2" x 2 3/8"
(349 mm x 419 mm x 60 mm)
Dent and rust proof platform

Product Footprint

14 3/8" x 21 1/8" x 53 5/8"
(359 mm x 537 mm x 1362 mm)

Product Weight

22.4 lb / 10.2 kg

Digital Height Range

3ft 11.2" – 6ft 7.5" / 120 cm – 202 cm
Resolution: 1 mm / 0.1"

Mechanical Height Range

11 3/4" – 79 1/2" / 30 cm – 202 cm
Graduation: 1/16" / 1 mm

Portability

2 wheels, wide 8" wheelbase for ease of movement

Warranty

2-year limited manufacturer warranty

ORDERING INFORMATION

MODEL	DESCRIPTION
502KL	Eye-Level digital scale with digital height rod, LB & KG
502KG	Eye-Level digital scale with digital height rod, KG only
502KL-BT	502KL with built-in Pelstar® wireless technology
502KG-BT	502KG with built-in Pelstar® wireless technology
500HB	Optional handlebar accessory for 502 series of scales
SS-502KL	ScaleSurance® 2 year extended warranty for 502 series*
ADPT31	Optional adapter for 502 series of scales
PLUGSET30	Plug adapters for use outside North America
BTKIT-01	Wireless upgrade kit for use with Welch Allyn® monitors
BTKIT-03	Wireless upgrade kit for use with pre-approved devices and Windows®/ Linux® PCs, USB dongle included
BTKIT-100	Wireless upgrade kit for use with pre-approved devices and Windows®/ Linux® applications, USB dongle not included
BTPRINTERKIT-01	Wireless thermal printer kit



Health o meter®
Professional



*Not available in all countries.