

**New!**  
Handheld model  
now available!

**Waldmann** **W**  
ENGINEERS OF LIGHT

---

## OPTICLUX HANDHELD

### CORDFREE LED MAGNIFICATION



#### Organic design. Outstanding quality of light.

The new OPTICLUX Hand magnifier luminaire is soon to be one of the medical practitioner's most valuable devices. Designed with first-rate performance and ergonomic comfort in mind, OPTICLUX Hand offers excellent color rendering properties, a large daylight quality light field, as well as

stepless dimming at your fingertips. In addition, the powerful Wood's light feature is ideal for veterinary or any dermatological task. Controls are easily accessible and a long battery life makes this a maintenance-free, dependable tool of the trade.



**Ergonomic**  
Optimal handling for right- and left-handers



**Easy to handle**  
Keypad on the handle with stepless dimming controls



**Outstanding light**  
Unique light output with 8 000 lx at 0.15 m for a true-to-life color rendering

The new OPTICLUX hand-held magnifier is the perfect choice for challenging visual tasks that require maximum accessibility. This is thanks to its large, dimmable daylight quality, its excellent color rendering characteristics and its long battery life. Medical professionals will be excited to work with the powerful Wood's Light in addition to daylight mode.



**Battery powered**  
Up to 3 hours runtime, convenient charging through USB-C port



**Optional charging base**  
Robust design with high quality baseplate. USB integrated charging cable included.



**Optimally protected**  
High quality felt cover protects the lens from scratches and sun exposure





#### OPTICLUX and Wood light.

The special mode for dermatologists, veterinarians and forensic medicine combines two lighting modes in one device: daylight white illumination and Wood light. With the flip of a switch, the practitioner is able to use the right tool for whatever task needs to be done. A variety of examination types can be performed with just one light. Other Wood luminaires require the room to be completely darkened. Thanks to powerful BLB LEDs, the Opticlux does not require this.

#### OPTICLUX HAND at a glance

- Color temperature daylight white 6,500 K
- Acrylic glass lens of 3.5 dioptries
- Color rendering Ra = 95
- Homogeneous, bright field of light \* (8,000 lx / 0.15 m)
- Continuously dimmable
- With Wood light: Maximum value at 365 nm
- Flexible charging via USB interface
- > 3 hours battery life
- Optional handheld model USB charging station
- Maintenance-free
- Ergonomic: perfect handling for right- and left-handers
- Optimally protected: with a high-quality felt case
- Tested medical device EN 60601-1
- 5-year warranty on all electronic components

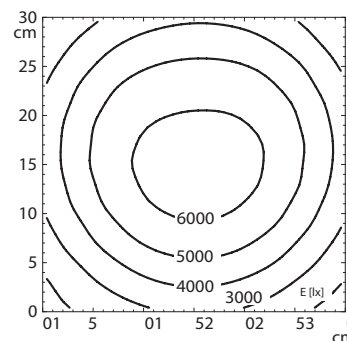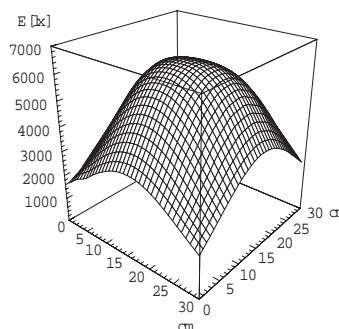

#### ORDERING INFORMATION

Model	Mounting Type	Special Features	Part No
Opticlux LED 10-2, UV - Handheld	Handheld	Wood light	D16067000
Opticlux USB charging station - accessory	Table or wall	USB integrated cable included	D16173000

**Waldmann** **W**

ENGINEER OF LIGHT.