


# TotalPure™ Centrifuge Tubes

## Metal-Free 15mL and 50mL Conical Tubes


- Ultra-pure polypropylene with less than 1ppm metal content
- 20,000×g maximum centrifugation rate
- Hi-Viz™ graduations with maximum freeze capacity line
- Resealable zip bags of 25 tubes, or individually wrapped
- E-Beam sterilized (SAL 10<sup>-6</sup>)

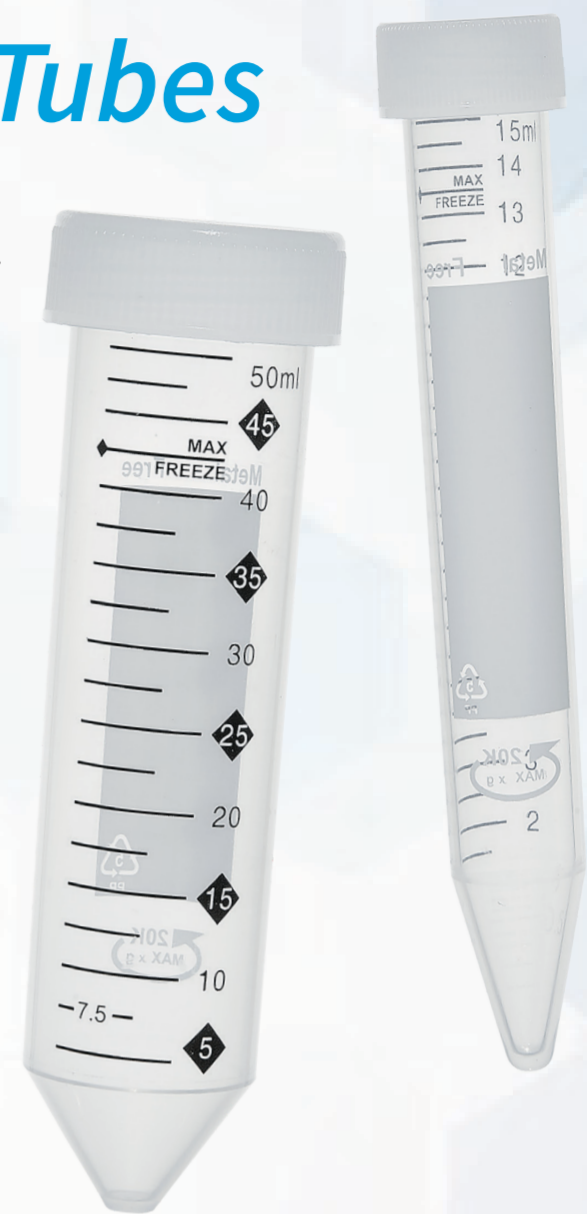

TotalPure™ Centrifuge Tubes are made with a proprietary purified polypropylene blend that reduces the metallic content present in most plastic blends to less than 1ppm. They also feature an array of premium performance and visual features.



The 20,000×g centrifugation rating is marked clearly on the tube. Hi-Viz™ graduations are easily identified at any angle, and feature a distinct indicator for the maximum fill level for freezing solutions. These tubes are recyclable with recycling information printed on the corner of the tubes' white marking area. HDPE caps are easy to open one-handed without compromising on the 95kPA leak-proof seal.

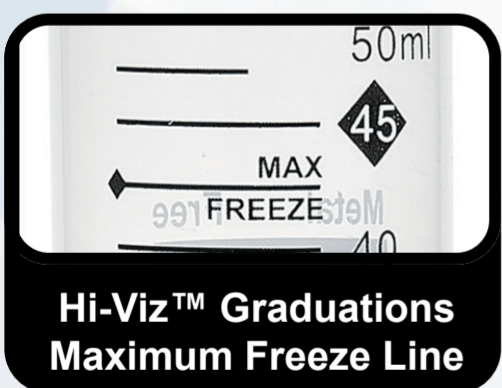

Tested to be certified free of: calcium, magnesium, zinc, iron, manganese, copper, aluminum, silicone, nickel, vanadium, sodium, phosphorous, cobalt, chromium, potassium, lithium, lead, selenium, cadmium, mercury, and arsenic.



CAT. NO.	DESCRIPTION	QTY
C2901	15mL Metal Free, Hi-Viz Graduations, Sterile, 20 resealable bags of 25	500 / case
C2915-W	15mL Metal Free, Hi-Viz Graduations, Sterile, Individually Wrapped	250 / case
C2903	50mL Metal Free, Hi-Viz Graduations, Sterile, 20 resealable bags of 25	500 / case
C2950-W	50mL Metal Free, Hi-Viz Graduations, Sterile, Individually Wrapped	250 / case