

## Power-Glide, Hi-Lo, Treatment Table with Adjustable Backrest, Footrest & Stirrups

### MODEL 87360

- One of the most affordable power tables in its class
- Features an adjustable backrest, friction-lock, footrest and stirrups
- Pneumatic backrest is adjustable from either side and raises to a maximum of approximately 80°
- Footrest is adjustable from either side and raises to approximately 75° from flat
- Retractable, stirrups swing out to a maximum of 35°
- Smooth-glide, electric power height adjustment
- Ultra heavy-duty base with solid, one-piece, steel ends
- Easy-clean, ABS plastic, center shroud hides lift components and provides a clean look
- Rigid, one-piece, laser-cut steel top frame
- Durable gray powder-coated steel frame
- Self-lubricating pivot points
- Smooth 110 volt electric height adjustment
- Power supply and motor meets UL 60601-1 and CAN/CSA C22.2 601.1
- Fast easy motor access via ABS side panels
- Hands-free low voltage foot control
- Seamless, rounded corner, upholstered top
- 2" ultra-firm padding
- Low height may help achieve ADA tax credit
- Customize with options
- Made in the USA by Clinton Industries
- 500 lbs. (227 kg) load capacity under normal use
- 3-year limited factory warranty from Clinton Industries

*power-glide™*  
Hi-Lo Treatment Tables by Clinton



California residents please be advised of the following, pursuant to Proposition 65:

**1**  
**WARNING:** This product can expose you to chemicals, including chromium, which is known to the State of California to cause cancer or reproductive toxicity. For more information go to [www.P65Warnings.ca.gov](http://www.P65Warnings.ca.gov).

### 87330 SPECIFICATIONS

Overall Size	76" L x 27" W x 18.5–34" H
Backrest	38" L x 27" W
Footrest	13" L x 18" W

### Upholstery Colors



## OPTIONS

### Positioning Straps Set

- 2" x 92" strap is secured with quick hook and loop Velcro® fasteners
- 2 straps per set (92" x 2")

080



### Safety Rails (pair)

- 3 height positions
- Plunger height adjustment

082



### Locking Casters

- Four, heavy duty, rubber wheel, easy-roll, locking casters

087



### 030 Paper Dispenser


- Accommodates paper rolls from 18" (45.72cm) to 21" (53.34cm) wide

030

### Paper Cutter

- Clear plastic with elastic strap and buckle

040

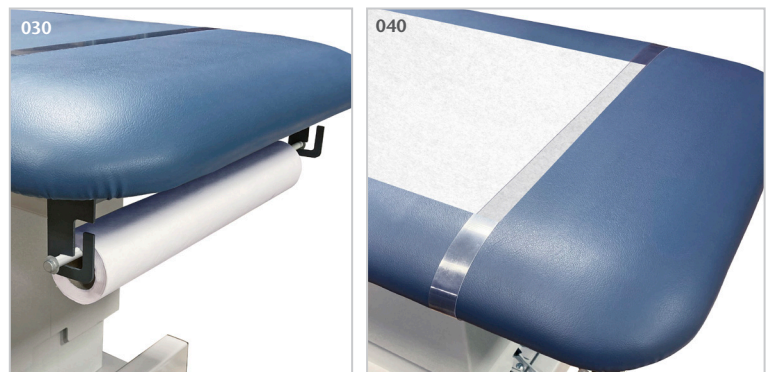

### Paper Dispenser/Cutter Combo

- Paper dispenser and paper holder ordered as one option

034