



space station

two tier back table

CDS-3060 Dimensions
 Main Table 30"D x 60"W x 34"H
 Top Shelf 20"D x 60"W



CDS-3072 Dimensions
 Main Table 30"D x 72"W x 34"H
 Top Shelf 20"D x 72"W


A specialty back table designed with additional room for large cases such as orthopedics, spinal fusions, neuro and craniotomies. The tuck-away cantilevered shelf holds multiple trays and is angled for clear observation of instruments.

The Space Station is easily rolled into position using 4" diameter, heavy-duty ball bearing brake/swivel casters. To lock swivel caster, simply engage the single-pedal brake.


One Person Operation



Pneumatic assist enables one-person operation in raising or lowering the cantilevered shelf.

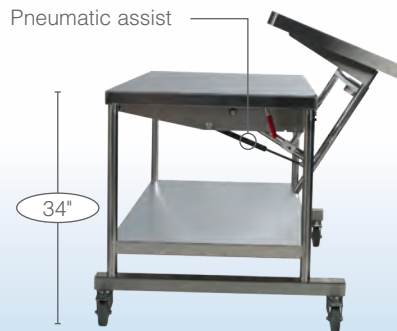

Shelf is easily raised when extra space is required, or lowered into storage position if not needed.

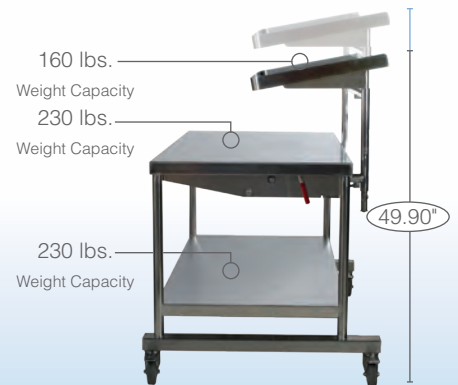

Cantilevered shelf features a 10° slope and has an adjustable height from 49.9" to 52.9".

trusted fit **drapes**
 two-tier

Disposable Table Drapes

Pedigo Trusted Fit drapes feature heavy-duty three layer construction. 12 Drapes per case.

CDS-3060-DD 5' table (1.5m)

CDS-3072-DD 6' table (1.8m)

100% u.s.a. supply chain

