

STAINLESS STEEL 10 SCOPE CAPACITY CABINET, KEY LOCK



SureDry
SCOPE CABINETS with

Dri-Scope Aid®
TRACKER

Scopes not included

SCSS8030TD-T2310D

FEATURES

- Positive pressure fully redundant fan system
- Unit features two intake fans and one exhaust fan
- Fan run time is fully adjustable and activated each time cabinet is closed
- HEPA filter removes 99.7% of particulates, .03 Microns
- Holds ten scopes with secure adjustable mounting system
- Replaceable HEPA Filter system without removal of scopes
- Medical grade electrical system
- Tempered glass double doors
- Multiple locking options available including badge entry
- Easy to clean removable drip tray
- Cabinets are delivered fully assembled
- Harloff has Industry Best 12 year limited manufacturers warranty for Sales to US and Canada

CABINET DIMENSIONS

- 93" H x 30" W x 24" D



*Also available in a
Powder Coat Finish*

Includes Dri-Scope Aid®, industry leading internal lumen drying system

Please refer to Dri-Scope Aid® specification sheet for more information

HARLOFF
INNOVATING STORAGE FOR CARE

Engineered &
Manufactured
in the U.S.



The ultimate in Endoscope drying systems is now available! SureDry Scope Cabinets feature an external HEPA filtered positive pressure drying system. Dri-Scope Aid® is the industry leading internal lumen drying system. This combination provides facilities with the assurance that they will be able to meet all current and anticipated future guidelines for external and internal drying of endoscopes, providing their patients the best possible assurance of infection prevention.

DRI-SCOPE INTERNAL LUMEN DRYING SYSTEM W/TRACKER




SureDry
SCOPE CABINETS with

Dri-Scope Aid®
TRACKER

The Dri-Scope Aid® Cabinet supplies HEPA-filtered air to all internal channels to aid in drying endoscopes with an automated drying cycle for improved turn-around time, and a continuous venting cycle for overnight storage. The Dri-Scope Aid® system in the SureDry Scope Cabinets provides a complete drying system for your facility.

TRACKER FEATURES & BENEFITS

- Digital display maintains a list of scopes by drying status, with displays of scope processing times
- Intelligent automated scheduler continuously processes scopes as long as they remain checked into the system
- Includes default or customizable settings and drying options unique to each endoscope
- Scope tracker maintains records of scope drying events, user access, check-ins and check-outs
- Color touch screen access by authorized personnel is completed by passcode, with adjustable session timeout settings and options for passcode retrieval if forgotten
- Various adjustable timing options controlled through the color touchscreen including drying, pre-drying and venting time settings, hang time limits, endoscope removal/load delay timing, etc
- At a glance, display shows scope inventory and reprocessing status including the time until scope expiration requiring further reprocessing
- Slot display order settings displayed on the status screen per user preferences
- Configurable replacement times for tubing and pump filter, with overdue warnings/indicator lights
- Downloadable export files available for endoscope management, user lists and more

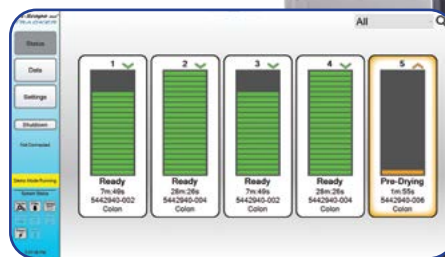

HARLOFF
INNOVATING STORAGE FOR CARE

Engineered &
Manufactured
in the U.S.



**Optional dimmable LED
lighting also available!**

10 Reasons to Choose **HARLOFF** for Your Scope Cabinets

A Smart Choice to Meet & Exceed Current SGNA, AORN and AAMI Guidelines



1

Industry Best Positive Pressure HEPA Filtered Drying System. Designed to add Dri-Scope Aid® now or later with easy “Plug & Play” installation.

Proudly Made in the USA

6



Harloff products are proudly **Designed, Engineered and Manufactured in the United States** at our Colorado Springs, CO facility.

LED Lighting Option

2




Brighten your cabinet and more easily inspect scopes with fully dimmable optional overhead LED lighting which triggers on/off when you open and close the door.

Easy Access to HEPA Filters

7

Change HEPA filters **without opening the cabinet.** Save time and ensure no cycle interruptions.

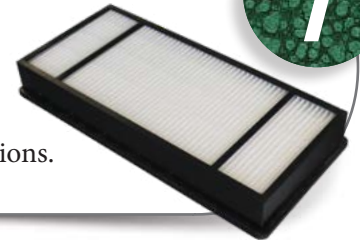

Keyed or Keyless Locks

8

Options include key locks, basic electronic push button locks, advanced keypads and proximity readers. View locking options online at www.harloff.com/locking-options/



Disposable Connectors

3

Save time and ensure compliance! Use single-use connectors for all internal channels of the endoscope. Part of the Dri-Scope Aid® package.



Hooks and Holders for Specialized Scopes

4

Major brands to include **Fuji, Olympus, Pentax, Karl Storz and others.**



Finish Options and Colors

9

Stainless Steel, Heat Sealed Powder Coat (High Gloss or Textured Low Gloss inside & out, view online at www.harloff.com/finish-options/) or upgrade to an **antimicrobial finish.**

12 Year Warranty

5

Our Manufacturer's Warranty is Unmatched in the industry. We build our cabinets to last and we stand behind them.

Mobile Scope Cabinet Kit

10

Convert any size of scope cabinet into a mobile version. Includes 4 3" casters and steel mounting plate.

