

DASH COAG

SET-AND-LOCK STAT COAG CENTRIFUGE

Part # 00-676-009-001



Reduce turnaround time (TAT) and produce high-quality coags in minutes with Drucker's DASH Coag centrifuge. The user-friendly set-and-lock controls streamline STAT lab operations and eliminate errors. Simply select your desired runtime from the DASH Coag's three preset cycles and spin with the press of a button.

Easy to Monitor

Prevent forgotten samples & reduce TAT with DASH Coag's LED lid lighting system. Bright lights indicate when the centrifuge is ready to load (off), running (on), or done (flashing) and continue flashing until the lid is opened.

Easy to Use

Stop wasting time adjusting settings or navigating complicated menus. The DASH Coag's arrives pre-programmed with the most common STAT settings and can spin tubes up to 16 x 100 mm or 10 mL with the same rotor and carriers.

Eliminate Wait Time

Further shorten TAT by using a centrifuge array in accessioning to eliminate queuing. The array takes advantage of the DASH Coag's compact size to allow samples to be spun as soon as they arrive. Its light weight and small footprint allow it to be easily moved throughout the lab to address bottlenecks.

Safe and Durable

Enjoy years of worry-free operation with the brushless motor and rugged, reinforced carbon fiber components. Enhance lab safety with imbalance protection and a clear, shatterproof lid, which locks until the centrifuge has stopped completely.


A Great Value

Each centrifuge comes ready to run with fixed angle rotor and 2-year warranty included in the affordable price.


Fixed Angle Rotor

The rotor remains at a 45° angle during processing for optimal PPP formation.

Loading Angle

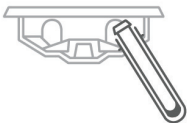

Rotation Angle

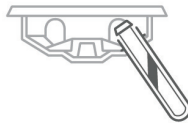

<p>Applications</p> <p>STAT Coag, PPP</p>	<p>Max Speed</p> <p>6600 Max RPM 4400 xg Max RCF</p>	<p>Max Capacity</p> <p>12 x up to 16 x 100 mm or 10 mL</p>	<p>Included</p> <p>Fixed Angle Rotor</p>	<p>Warranty</p> <p>2 years Extended warranty options available</p>
--	---	---	---	---

DASH COAG



Fixed angle centrifugation for optimal coag production

Lid lights indicate cycle status: ready, running, done

3 default presets for common coagulation cycles

QR code for easy access to support videos, manuals, and more

Press to select cycle

Use presets or program your validated cycles using the 2-button interface

May be locked to a single cycle or customized at the lab

Free 30-Day Evaluation (U.S. Only)

Standard Preset Cycles	Time	RPM	G-Force
Cycle 1	2 min	6600	4400
Cycle 2	3 min	6600	4400
Cycle 3	5 min	5200	2700

Part Numbers	
Centrifuge	00-676-009-001
1-1.5 mL Microtainer Tube Adapter (6 qty)	03-1-0007-0059K
1 yr Warranty Extension	03-0-0007-0003
2 yr Warranty Extension	03-0-0007-0004
3 yr Warranty Extension	03-0-0007-0005

(adapters available for a variety of tube sizes)

Specifications	
Sound Level	61 dBA
Max RPM	6600 RPM
Max G-Force	4400 xg
Dimensions (W x D x H)	12 x 14 x 9 in 30.5 x 35.6 x 22.9 cm
Weight	12 lbs (5.4 kg)
Power	100-240 (+/- 10%) - 50/60 Hz

