


BREATHEDRY™

BreatheDRY™ BridgeAIR™ Absorbent Protection Pads


BreatheDRY™ Breathable Absorbent Pads are a white, soft flat sheet with an absorbent layer and a waterproof breathable underside that stops fluid from getting through to equipment below and wicks away moisture to protect the patient's skin—designed to be breathable and to manage climate control and preserve skin integrity. The BreatheDRY™ Absorbent Protection Pads provides a breathable, latex free absorbent layer that can be placed on top of the BridgeAIR™ mattresses to stop fluids from getting through and protect a patient's skin while providing comfort.

3 Layers of Protection:

1. Soft, non-woven air-permeable top layer is gentle on the skin of the very young to the very old, as well as vulnerable and fragile skin during periods of illness.
2. The super absorbent inner layer allows body fluids to be taken away from the skin. The fluids are locked in the inner layer to keep skin dry.
3. The breathable, waterproof backing assists with the containment of fluid for safe disposal.

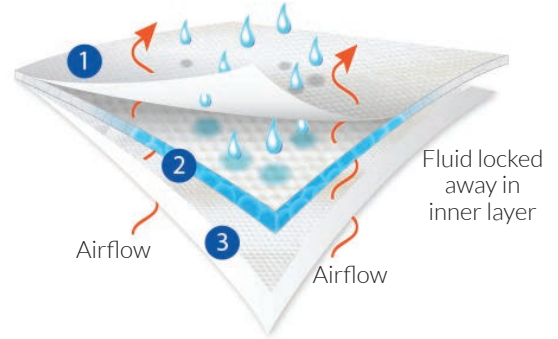


**99.87%
MOISTURE
IS LOCKED
IN AFTER 5
MINUTES!**

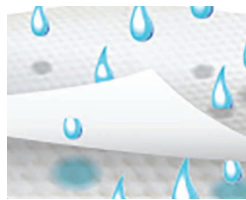

Moisture Management

Preserving the micro-climate, the BreatheDRY™ draws in moisture, locks it away from the skin and feels dry to the touch in just minutes.

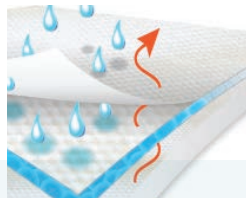

Breathable Backing

Permits air flow with air permeable backsheet. Compatible with all pressure redistribution surfaces, including low air loss mattress therapy. 20 times more breathable than traditional blue underpads.



High Level Fluid Handling Capacity

Absorbs 3.2 liters per m² (AWTA 62) which assists the containment and safe disposal through control of fluids.



Reduce Cleaning Costs

By protecting your bedding, airassisted patient transfer mats, or tables, equipment, and patient slings.

SIZES AND SPECIFICATIONS

BREATHEDRY™ PROTECT FULL, 39" X 78", BRIDGEAIR™ ABSORBENT PROTECTION PAD

Absorbent Area: 35" x 70"
Maximum Absorbency: 6.4 liters
Product Weight: 8 oz.

SKU: BH-BD8040
QTY: Individually sealed, 20/cs

BREATHEDRY™ STANDARD 31" X 35", BRIDGEAIR™ ABSORBENT PROTECTION PAD

Absorbent Area: 27" x 31"
Maximum Absorbency: 2.3 liters
Product Weight: 4 oz.

SKU: BH-BD3135
QTY: Individually sealed, 60/cs

BREATHEDRY™ STANDARD 37" X 56", BRIDGEAIR™ ABSORBENT PROTECTION PAD

Absorbent Area: 31" x 47"
Maximum Absorbency: 3.07 liters
Product Weight: 4 oz.

SKU: BH-BD3756
QTY: Individually sealed, 60/cs

BREATHEDRY™ PLUS 48" X 39", BRIDGEAIR™ ABSORBENT PROTECTION PAD

Absorbent Area: 39" x 33"
Maximum Absorbency: 3.8 liters
Product Weight: 4 oz.

SKU: BH-BD4839
QTY: Individually sealed, 30/cs

