



TRUEBOND

MULTIPURPOSE VETERINARY CENTRIFUGE

Part # 00-376-009-005

Using the TrueBond multipurpose veterinary centrifuge, vet laboratories can now spin microhematocrit tubes, fecal, blood, urine, and more in one centrifuge.

Spin Any Tube

TrueBond is designed to spin all your tubes, from microhematocrit to urine and everything in between. Includes all necessary adapters.

Easy to Use

Minimize training time with the user-friendly interface. The fully digital display makes it easy to customize cycle parameters on the fly, or set and select one of the 10 programmable cycles.

Easy to Monitor

Process samples more efficiently with TrueBond's LED lid lighting. Once the unit is started, the display shows remaining spin time.

Safe and Durable

Enjoy years of safe, worry-free operation with the brushless motor, imbalance protection, and a clear, shatterproof lid, which locks until the centrifuge has stopped completely.


A Great Value

Each TrueBond centrifuge is ready to run with a fixed angle rotor, two sets of carriers, microtainer adapters, microhematocrit adapters and reader card, and a 2-year warranty included in the affordable price.



Applications	Max Speed	Max Capacity	Included	Warranty
Blood, Urine, Fecal, Microhematocrit	3900 Max RPM 1850 xg Max RCF	6 x blood, urine, and fecal tubes Up to 16 x 125 mm or 15 mL 8 x Microhematocrits	Fixed Angle Rotor 6 x 16 x 100 mm Blue Carrier 6 x 16 x 125 mm Orange Carrier 6 x 1-1.5 mL Microtainer Adapter 2 x Microhematocrit Adapters Veterinary Microhematocrit Reader Card	2 years Extended warranty options available

TrueBond

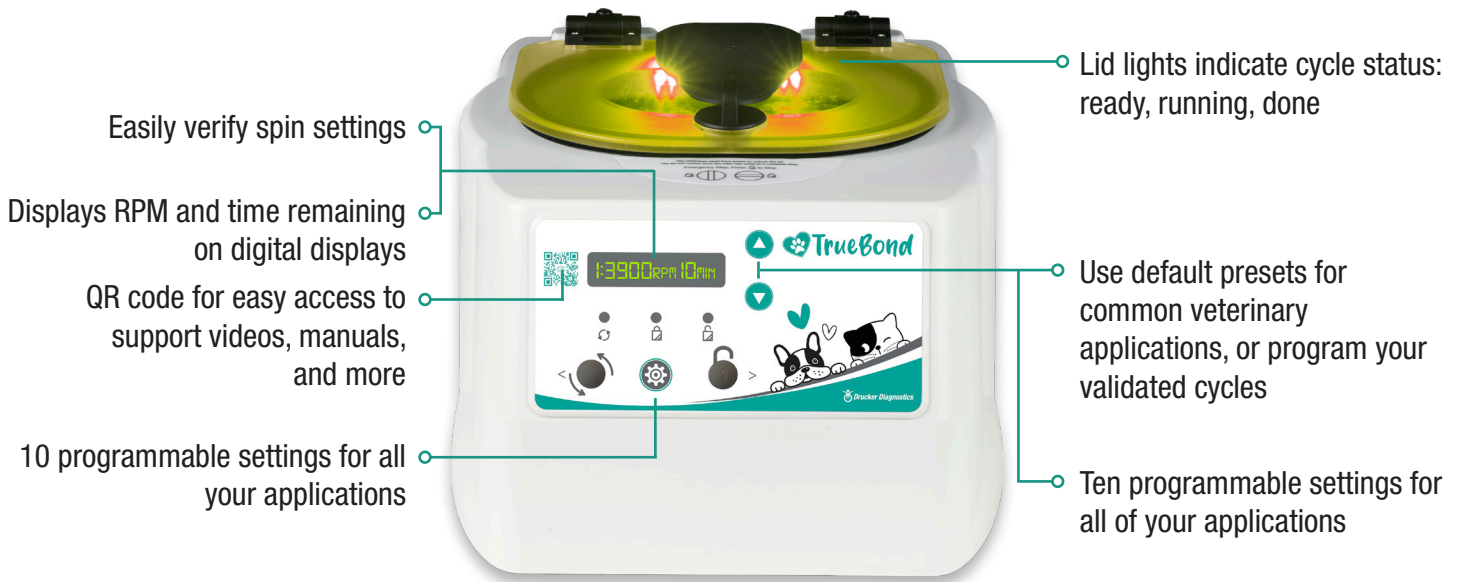

Free 30-Day Evaluation (U.S. Only)



Standard Preset Cycles	Time	RPM	With Blue Carriers	With Orange Carriers
			(Up to 16 x 100mm or 10 mL tubes)	(Up to 16 x 125 mm or 15 mL tubes)
Cycle 1 (Chemistry)	10 min	3900	1650 xg	N/A
Cycle 2 (Urine Sediment)	5 min	2000	N/A	500 xg
Cycle 3 (Veterinary Fecal)	5 min	1600	300 xg	N/A
Cycle 4 (Microhematocrit)	6 min	3800	1550 xg	N/A

Part Numbers

Centrifuge	00-376-009-005	
Included	16 x 100 mm Carrier, Blue (6 qty)	7713079K
	16 x 125 mm Carrier, Orange (6 qty)	7713044K
	1-1.5 mL Microtainer Tube Adapter (6 qty)	03-1-0007-0059K
	Microhematocrit Adapter (2 qty)	02-004-0-0027K
	Veterinary Microhematocrit Reader Card	02-0006-0-0046
Optional	1 Year Warranty Extension	03-0-0007-0001
	2 Year Warranty Extension	03-0-0007-0002
	3 Year Warranty Extension	03-0-0007-0006

(adapters available for a variety of tube sizes)

Specifications

Sound Level	64 dBA
Max RPM	3900 RPM
Max G-Force (Blue Carriers)	1600 xg
Max G-Force (Orange Carriers)	1850 xg
Dimensions (W x D x H)	12 x 14 x 9 in 30.5 x 35.6 x 22.9 cm
Weight	12 lbs (5.4 kg)
Power	95 to 253 V – 50/60 Hz

