


AlarmKIT

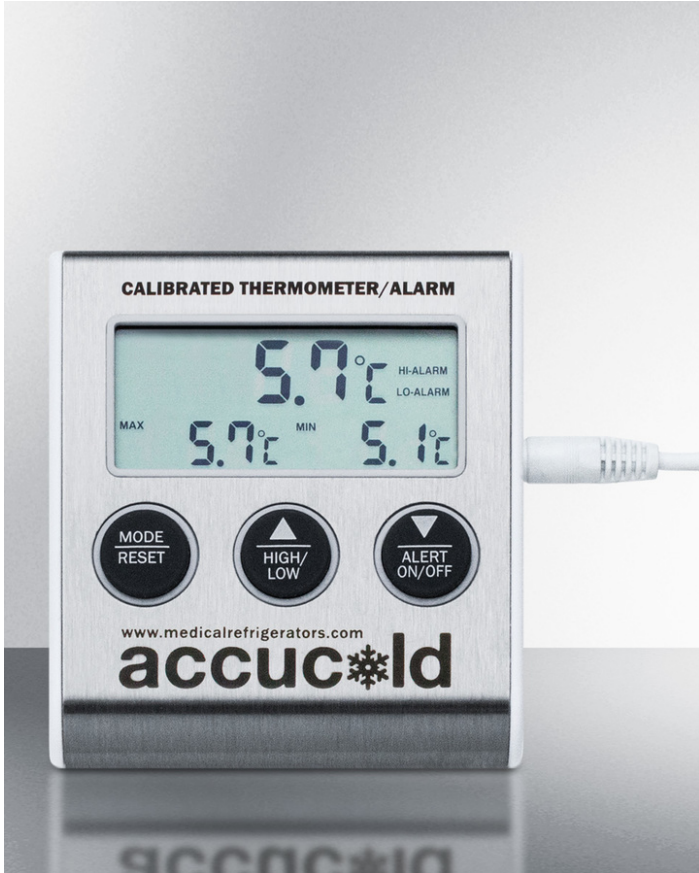

3" x 2.75" x 0.5" (H x W x D)

High/low temperature alarm with NIST calibrated temperature readout

Highlights:

Temperature probe is calibrated to NIST standards in our ISO/IEC 17025:2017 accredited laboratory

Audible and visual alarms alert users if temperature rises or falls below the specified approved range

LCD display lets you view the current, minimum, and maximum temperature to the nearest tenth of a degree

Product Features:

NIST calibrated	Temperature probe is calibrated in our ISO/IEC 17025:2017 accredited laboratory
Includes NIST certificate of calibration	All alarms ship with NIST calibration certificate to save for your facility's records
Audible/visual alarm	Alarm sounds and screen flashes if the unit temperature rises above or falls below the high and low set points
Protected temperature sensor	Temperature sensor is enclosed inside glycol-filled bottle to protect it from being affected by rapid temperature changes that may occur when the unit's door is opened
Reliable temperature reporting	Accurate within .5° C (1° C when installed on units operating below -30°C or above 70°C; outside specified range, accuracy is 1°C)
Temperature display	External digital display shows the current, minimum, and maximum temperatures to the nearest tenth of a degree
Viewable in Celsius or Fahrenheit	Rear slide switch lets you switch the display from Celsius to Fahrenheit
Battery-powered	Includes a replaceable AAA lithium battery
Low battery visual alert	Battery symbol flashes on the LCD screen when battery needs to be replaced

AlarmKIT Specifications:

Overview	
Nominal Height	3.0" (8 cm)
Height of Cabinet	3.0" (8 cm)
Nominal Width	3.0" (8 cm)
Width	2.75" (7 cm)
Nominal Depth	1.0" (3 cm)
Depth	0.5" (1 cm)
Cabinet	Stainless Steel
Weight	1.0 lbs. (0 kg)
Parts & Labor Warranty	90 Days
UPC	761101016184

