

AFS23ML



83.75" x 27.5" x 31" (H x W x D)

Performance series Pharma-Lab 23 cu.ft. all-freezer in stainless steel

Highlights:

Optimized forced air cooling for excellent temperature stability and uniformity with rapid recovery after door openings and product loading

Audible/visual alarms for high/low temperature excursions

Temperature viewable in Celsius or Fahrenheit to the nearest degree

Product Features:


Purpose-built design	Designed and engineered to store vaccines and other sensitive materials under stable temperature conditions
Adjustable temperature	Interior is pre-set to -25°C, with an adjustable control range from -30 to -15°C
Temperature alarm	Visual and audible alarms alert if the interior temperature rises or falls out of set range
Complete stainless steel construction	Interior and exterior are entirely constructed from stainless steel for lasting durability and easy cleanup
Forced air cooling	Top-mounted commercial fan ensures rapid temperature recovery
Automatic defrost with adjustable cycle defrost	Adjustable cycle defrost operation ensures the unit avoid temperature spikes and ice buildup
Self-closing door	Closes to save energy but will remain open when the door is opened at a 90° angle for easier loading
Open door alarm	Alarm will beep 3 times if the door is left open for 30 seconds, 5 times at 60 seconds, and will sound continuously if the door is open for more than 5 minutes
Commercial grade casters	Ships with a set of casters (2 locking) for easy mobility as needed




AFS23ML Specifications:

Overview	
Nominal Height	84.0" (213 cm)
Height of Cabinet	83.75" (213 cm)
Nominal Width	28.0" (71 cm)
Width	27.5" (70 cm)
Nominal Depth	31.0" (79 cm)
Depth	31.0" (79 cm)
Depth with Handle	32.75" (83 cm)
Depth with door at 90°	56.5" (144 cm)
Capacity	23.0 cu.ft. (651 L)
Defrost Type	Frost-Free
Door	Stainless Steel
Cabinet	Stainless Steel
US Electrical Safety	UL
Canadian Electrical Safety	ULC
Sanitation	ETL-S
Amps	8.0
Voltage/Frequency	115 V AC/60 Hz
Weight	273.0 lbs. (124 kg)
Parts & Labor Warranty	2 Years
Compressor Warranty	5 Years
UPC	761101062310
Freezer	
Door Swing	RHD
Reversible	No
Interior Height	51.5" (131 cm)
Interior Width	23.5" (60 cm)
Interior Depth	27.0" (69 cm)
Shelf Type	Wire
Shelf Qty	4
Adjustable Shelves	Yes
Thermostat Type	Digital
Interior Light	Yes
Fan Type	Interior
Refrigerant Type	R290
Refrigerant Amount	3.5 oz.
Wheels	4
Temperature Range	-30 to -15°C

Sound Level

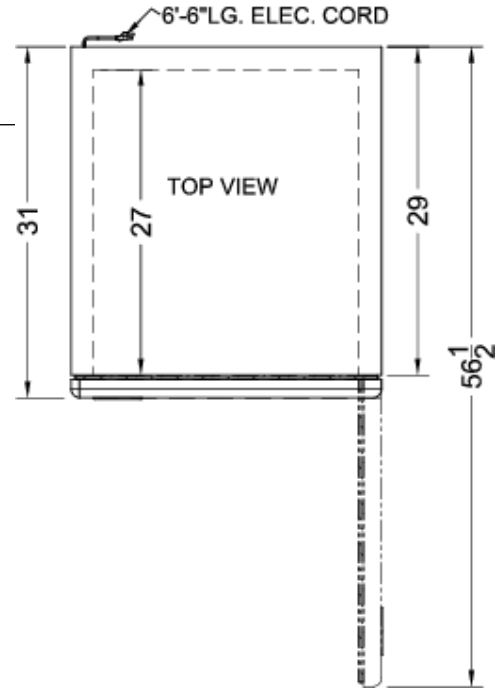

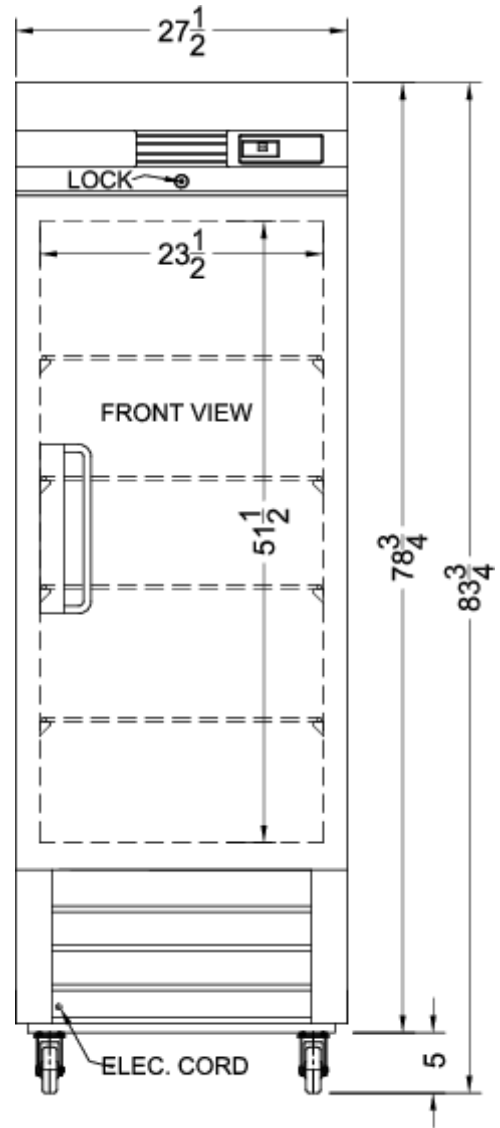

AFS23MLLH



83.75" x 27.5" x 31" (H x W x D)

Performance series Pharma-Lab 23 cu.ft. all-freezer in stainless steel with left hand door swing

Highlights:

Intelligent microprocessor digital temperature controller

Optimized forced air cooling for excellent temperature stability and uniformity with rapid recovery after door openings and product loading

Audible/visual alarms for high/low temperature excursions

Product Features:

Purpose-built design	Designed and engineered to store vaccines and other sensitive materials under stable temperature conditions
Adjustable temperature	Interior is pre-set to -25°C, with an adjustable control range from -30 to -15°C
Temperature alarm	Visual and audible alarms alert if the interior temperature rises or falls out of set range
Complete stainless steel construction	Interior and exterior are entirely constructed from stainless steel for lasting durability and easy cleanup
Forced air cooling	Top-mounted commercial fan ensures rapid temperature recovery
Automatic defrost with adjustable cycle defrost	Adjustable cycle defrost operation ensures the unit avoid temperature spikes and ice buildup
Self-closing door	Closes to save energy but will remain open when the door is opened at a 90° angle for easier loading
Open door alarm	Alarm will beep 3 times if the door is left open for 30 seconds, 5 times at 60 seconds, and will sound continuously if the door is open for more than 5 minutes
Commercial grade casters	Ships with a set of casters (2 locking) for easy mobility as needed



AFS23MLLH Specifications:

Overview	
Nominal Height	84.0" (213 cm)
Height of Cabinet	83.75" (213 cm)
Nominal Width	28.0" (71 cm)
Width	27.5" (70 cm)
Nominal Depth	31.0" (79 cm)
Depth	31.0" (79 cm)
Depth with Handle	32.75" (83 cm)
Depth with door at 90°	56.5" (144 cm)
Capacity	23.0 cu.ft. (651 L)
Defrost Type	Frost-Free
Door	Stainless Steel
Cabinet	Stainless Steel
US Electrical Safety	UL
Canadian Electrical Safety	ULC
Sanitation	ETL-S
Amps	8.0
Voltage/Frequency	115 V AC/60 Hz
Weight	273.0 lbs. (124 kg)
Parts & Labor Warranty	2 Years
Compressor Warranty	5 Years
UPC	761101066738
Freezer	
Door Swing	LHD
Reversible	No
Interior Height	51.5" (131 cm)
Interior Width	23.5" (60 cm)
Interior Depth	27.0" (69 cm)
Shelf Type	Wire
Shelf Qty	4
Adjustable Shelves	Yes
Thermostat Type	Digital
Interior Light	Yes
Fan Type	Interior
Refrigerant Type	R290
Refrigerant Amount	3.5 oz.
Wheels	4
Temperature Range	-30 to -15°C

Sound Level

