



## The Importance of Animal Scale Calibration

Rice Lake's certified calibration weights are ideal weights to use when calibrating animal scales. They ensure your scale is properly calibrated and ready for the requirements of weighing in a wide range of animal care settings

Rice Lake Certified Calibration Weights ensure our scale is properly calibrated and ready to weigh a wide range of animal needs.

### Weights for Animal Scale Calibration

Consumer weights such as barbells or similar items may not provide an accurate weight to calibrate your scale. For example, a 5-pound consumer weight, such as a barbell, may actually weigh between 4.75 and 5.25 pounds. Certified calibration weights have a manufactured and adjusted tolerance of +/-230 milligrams. Using a certified calibration weight will ensure your scale is calibrated correctly and presents accurate weighments for proper animal care.

Rice Lake Weighing Systems offers a wide range of certified calibration weights designed to verify the accuracy of your scales. Visit [www.ricelake.com/precision](http://www.ricelake.com/precision) to view our calibration weights

# Rice Lake VS-12 Dual Range



Model Rice Lake VS-12

Approvals



## Standard Features

- Large, easy-to-read, 0.75 in LCD display
- Dual-range capacity for increased precision from 0 to 10 lb
- Automatic precision change based on patient weight
- Unit of measure locking
- Automatic shut-off for long battery life
- Motiontrap™ movement compensation technology
- Built-in measuring tape 0 to 23.5 in × 0.125 in (0 to 60 cm × 0.25 cm)

## Specifications

### Capacity:

0 to 10 lb × 0.05 oz (0 to 5 kg × 1 g)  
10 to 33 lb × 0.1 oz (5 to 15 kg × 2 g)

### Power:

12 VDC, provided by six AA lithium batteries (included) or AC adapter (optional)

### Battery Type:

Six Li-FeS<sub>2</sub> batteries

### Operating Temperature:

50 °F to 95 °F (14 °C to 35 °C)

### Display:

5-digit LCD display, 0.75 in (1.9 mm) digit height

### Approvals:

c-Met E113986

### Warranty:

Two-year limited

## Part Number/Price

Part #	Model #	Capacity	Price
204398	Rice Lake VS-12	0 to 10 lb × 0.05 oz, 10 to 33 lb × 0.1 oz (0 to 5 kg × 1 g, 5 to 15 kg × 2 g) 120V/230V	Consult

## Shipping Weight


Part #	Model #	L (in)	W (in)	H (in)	Package Wgt. (lb)	Dim Wgt. (lb)	Shipping Wgt. (lb)
204398	Rice Lake VS-12	18	28	11	25	40	40

## Options/Accessories

Part #	Description	Price
165138	AC adapter, input voltage 100 to 240 V, c-Met, medical grade	Consult

## Dimensions

	Inches	Centimeters
A	24.0	60.8
B	19.0	48.2
C	6.8	17.1
D	17.0	43.2
E	3.5	8.9

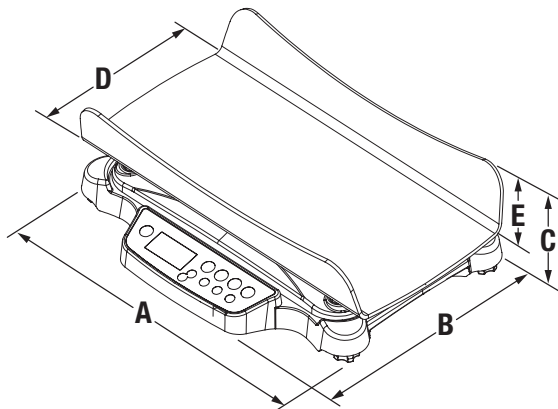

COMPANION  
ANIMAL SCALES