



C-Matt™ Prevention Plus Replacement Mattress

Top of the Line Mattresses at Manufacturer Direct Price

C-Matt Prevention Plus replacement mattresses are quality patient support surfaces that offer exceptional prevention of skin breakdown in addition to safety and comfort benefits. Our mattresses are custom designed to meet your every need — **patient comfort and safety, infection control, and easy cleaning and maintenance.** All C-Matt replacement mattresses are made in the USA by a Certified Woman Owned Small Business.



What C-Matt Prevention Plus Replacement Mattress Can Do for You

Our C-Matt Prevention Plus replacement mattress fits any standard flat deck bedframe. It is designed to prevent skin breakdown, effectively treat existing skin breakdown and improve patient comfort. Every detail works to enhance the safety and convenience of your patients, medical staff and environmental safety team.

Features & Benefits

Cover features include:

- **Patented seal.** Radio frequency (RF) welding technology and our patented “waterfall” zipper flap prevent fluids from penetrating the mattress.
- **Bi-directional stretch.** Breathable, two-way stretch fabric gives added comfort and durability. Also prevents “hammocking;” allowing for easier patient transfer.
- **Microclimate technology.** Polyurethane coated fabric absorbs body heat and perspiration and disperses it away from the body — keeping patients cool and dry, reducing the incidence of skin breakdown.
- **Full perimeter zipper.** A zipper that goes around the cover’s entire perimeter makes for easy removal for cleaning and maintenance.
- **Antimicrobial and hypoallergenic.** Our fabric is treated to inhibit bacterial growth and has proven to be hypoallergenic.
- **Medical-grade tri-laminate vinyl bottom.** The bottom of our cover is constructed of a polyester scrim surrounded by two layers of vinyl. This durable construction enhances the longevity of the mattress and helps it maintain its proper form. It is also antimicrobial, antifungal and fluid-proof.
- **Cover options.** All of our covers and liners are designed to help eliminate shear and friction; drastically reducing the risk for skin breakdown. Replacement covers are available.



Foam features include:

- **Multi-zone foam core.** High resiliency foam provides the right support in the right places. Body weight is redistributed to reduce pressure on high-risk body points. Patients immerse into the mattress and are comfortably supported while enveloped by its multi-zone foam. Meets Cal 117-2013 fire code.
- **Sloped convoluted heel section.** The generous sloped heel section reduces the need for adjunctive heel suspension. Extra soft convoluted foam in this typically high-pressure area provides additional comfort and protection from skin breakdown.
- **Firm perimeter side rails.** The perimeter of our mattress consists of heavy-duty foam that provides added security and support to patients. These supportive side rails provide a barrier to keep patients safely nestled toward the center of the mattress and prevent bottoming out when patients sit on the side of the mattress to get in and out of bed.
- **Shear liner.** Our foam is surrounded by a liner to eliminate friction between the foam core and the cover. The same microclimate technology used in our covers is also used in this liner to absorb and redistribute body heat and perspiration and to reduce shear and friction.
- **Fire barrier.** (Optional) Our barrier is in compliance with all state and federal fire safety standards. Meets Cal 129/Fed 1633.




Specialized foam zones for optimal pressure redistribution



Dedicated extra soft convoluted heel section



RF welded seams to prevent cross contamination

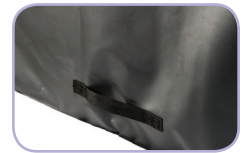

Top of the Line Mattresses at Manufacturer Direct Price

Specifications:

Specifications

Brand Name
Depth
Facility Type
Fire Barrier
Latex Free
Length
Mattress Core Material

Mattress Cover Material

Risk Level

Weight Capacity
Width

Values

C-Matt Prevention Plus
6"
Acute Care
Optional
Yes
80" and 84"
High Density & Visco Foam
Deep Cut Convolute Designed To Cool
Sloped Heel Section
B-Directional Stretch Urethane
Non-Slip Tri-Laminate Vinyl
High Risk Early Stage 3

500 lbs.
36"

Packaging Details

Unit of Measure

Each

Weight

35 lbs.

Package Dimensions

37 X 30 X 19