




# C-Matt™ StatiCare

*A High Performing Static Air Surface to Meet Your Facility's Skin Safety and Comfort Needs*

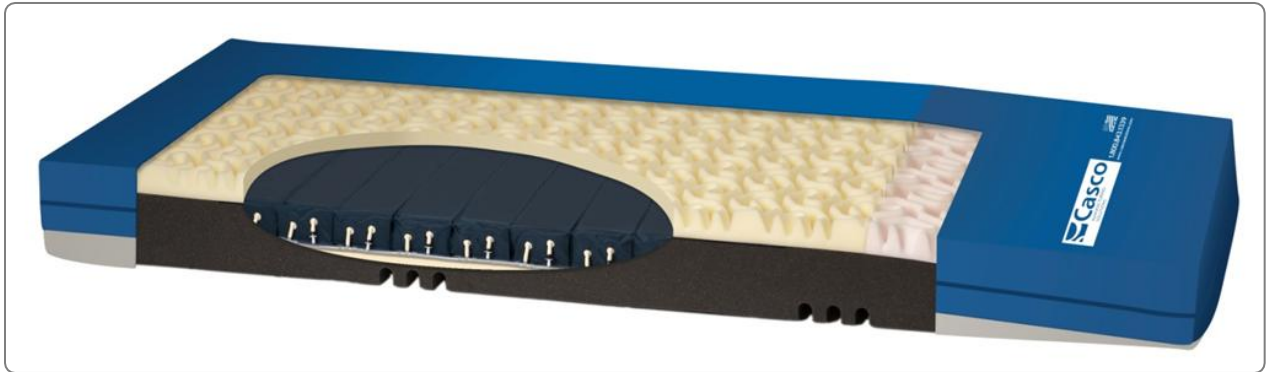

## A Unique Approach to Pressure Redistribution

### Best-in-Class Top Cover

*Breathable, Bi-directional Stretch Fabric  
Fully RF Welded to Prevent Fluid Ingress*

### Interconnected Air Cells for Even, Quiet, Non-Powered Pressure Redistribution

### Generous Dedicated Sloped Heel Section

### Unique Convoluted Foam Topper for Added Comfort

### Standard Additional Fire Barrier for Regulatory Compliance

### Patient Weight Limit: 500lb

## A Quality Replacement Surface System

### Interconnected Air Cell Section

- Eight horizontal, anatomically located air cells for maximum efficacy
- Requires no electricity, quiet
- Delivers dynamic response to patient movement

### Surface Foam Core

- High density, anti-microbial, medical grade foam
- Firm base and perimeter for patient safety
- Large dedicated heel section

### Replaceable Cover

- Anti-microbial and hypo-allergenic
- Easy to clean, stain resistant and durable
- Patented zipper-flap design to protect mattress core
- Full perimeter zipper for easy maintenance
- Replacement Covers available

### Standard (7 x 36 x 80 / 84) and Custom Sizes Available

### 5 year Warranty

Interconnected air chambers automatically adjust to movement



Dedicated extra soft convoluted heel section



RF welded seams to prevent fluid ingress



Meets Fire Safety Regulatory Requirements CAL117, CAL129

