

Surgical light

AIM-100



Limitless rotation, class-leading intensity and endless versatility.

Providing both amazing value and unparalleled aesthetics, Philips Burton's AIM-100 also boasts limitless rotation on ceiling mounts (a feature previously available only on OR type lights). Along with highest intensity illumination, AIM-100 delivers very smooth focus and superior maneuverability. A longer arm system ensures that either light can be positioned in the widest range of positions possible, making it ideal for examination applications.



AIM-100


- High-intensity of 63,000 lux (5852 fc) at 1 meter (per lighthead on double ceiling model)
- 3500 K color temperature offers natural and true coloring
- High CRI (Color Rendering Index) of 97
- Uses three quartz halogen bulbs (per lighthead on double ceiling model) – each rated with a 2000-hour bulb life; bulbs wired in parallel for continuous operation
- Compact 20" (51 cm) diameter lighthead
- Extremely smooth focus mechanism with center sterilizable handle
- 360° limitless arm-and-mounting-system-rotation around vertical axes (on single ceiling models)
- Cool operation for long hours of use
- 120V standard; 100 and 230 / 240V models available
- UL/IEC 60601-1, IEC 60601-1-2, EN/IEC 60601-2-41 and CAN/CSA C22.2 NO. 601.1 M90 certified
- CE marked
- 5-year limited warranty
- Assembled in USA

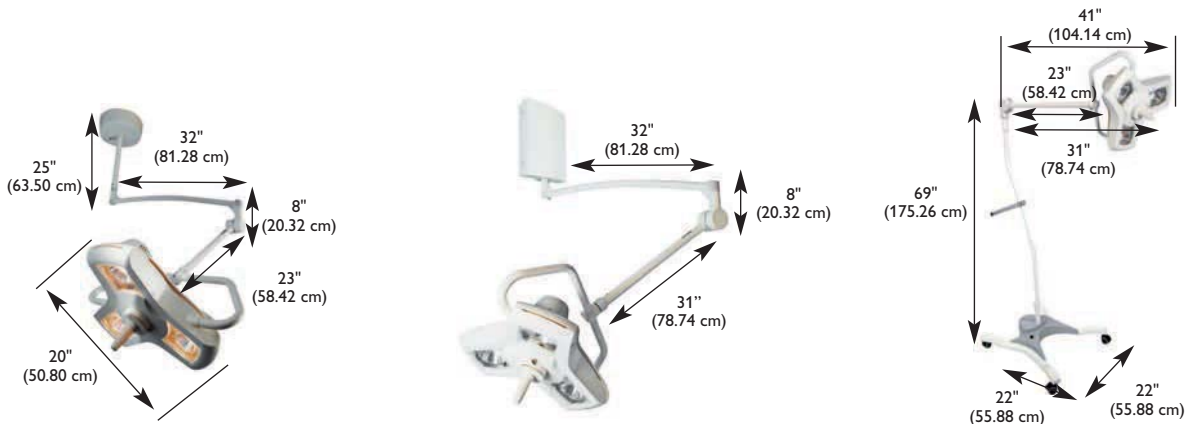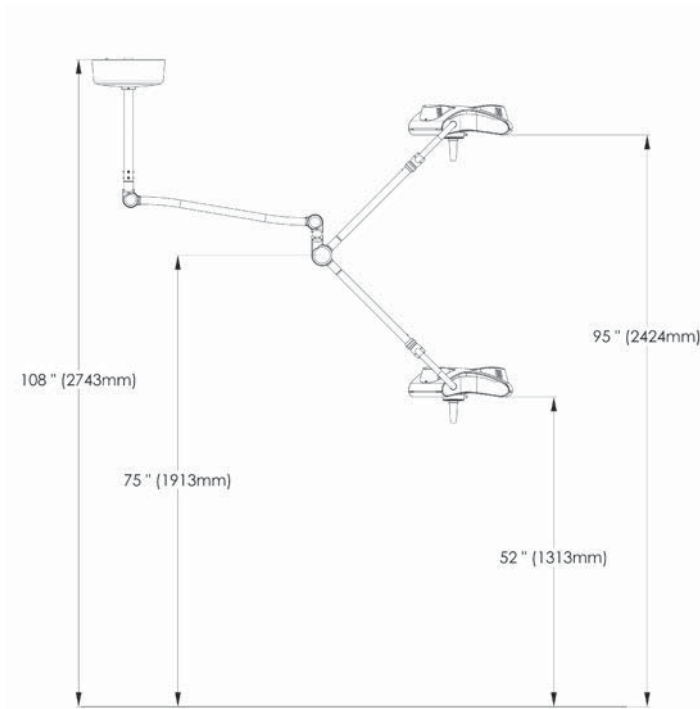



AIM-100 Double Ceiling

Philips Burton's AIM-100 Double Ceiling provides both amazing value and unparalleled aesthetics. Along with the highest intensity illumination, you can also expect a smooth focus and superior maneuverability. A longer arm system ensures that the light can be positioned in the widest range of positions possible, making the AIM-100 Double Ceiling ideal for surgery.

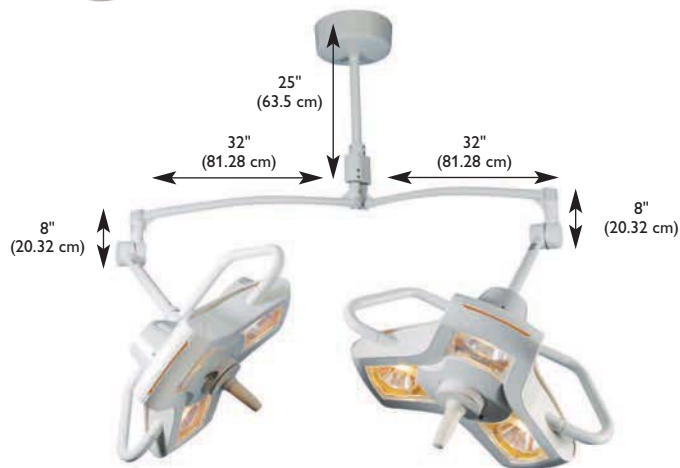

Ordering information

120V Model	230/240V Model	230V SP Model*	100V Model	Description
A100FL	A103FL	A102FL25	A101FL	Floorstand
A100SC	A103SC	A102SC25	A101SC	Single ceiling mount
A100DC	A103DC	A102DC25	A101DC	Double ceiling mount
A100W	A103W	A102W25	A101W	Wall mount
Accessories				
1017111PK				Set of 3 AIM-100 replacement bulbs
1017040				Sterilizable / autoclavable handle
0008100PK				Disposable handle covers, 25/box
1017084				42" drop tube

*SP Model refers to Spanish Language Manuals.