

# accucold®




## PHARMA-VAC PERFORMANCE

# DUAL COMPRESSOR COMBINATION REFRIGERATOR-FREEZER

Designed and purpose-built for pharmacy, medication, vaccination, and life science applications to support meeting CDC/VFC and Health Canada vaccine storage guidelines



- ✓ Independent cooling of each compartment
- ✓ Audible & visual high/low temperature alarms
- ✓ Optimized efficient performance for stable temperature
- ✓ Built-in readout for current & high/low temperature
- ✓ Anti-microbial handles to reduce germ spread
- ✓ Eco-friendly natural refrigerants for a greener footprint

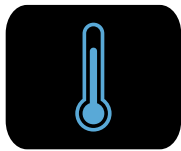

## ADVANCED TEMPERATURE CONTROL

- Intelligent microprocessor digital temperature controllers
- Improved temperature accuracy with buffered probe in glycol-filled bottles to simulate temperature of stored product
- Additional sensor reads interior air temperature for superior temperature control
- Adjustable operating control range: Refrigerator: from +2°C to +8°C Freezer: -30°C to -10°C
- Digital display of the min/max temperature (updated in 15 minute increments)
- Temperature viewable in Celsius and Fahrenheit to the nearest tenth of a degree
- Password protection of control parameters beyond setpoint
- Optimized forced air cooling for excellent temperature stability and uniformity with rapid recovery after door openings and product loading



## MULTIPLE ALARMS FOR SAFER MONITORING

- Audible/visual alarms for high/low temperature excursions
- Open door alarm sounds if the door is left ajar (default 4 minutes with selectable 1-99 minute range)
- Audible power failure and sensor failure alarms
- Mute function for alarms
- Remote alarm contacts in rear of the unit



## FULLY FEATURED INTERIOR FOR CONVENIENT STORAGE

- Energy efficient LED lighting with internal on/off switch
- Adjustable plastic-coated wire shelves for improved circulation
- Multiple shelf positions spaced at ½" to allow flexible adjusting to support smaller items
- Adjustable cycle defrost operation (Refrigerator) and static manual defrost (Freezer)
- ½" probe ports included to accept user-provided monitoring equipment



## QUALITY DESIGN FOR RELIABLE PERFORMANCE

- Hospital grade cord with a green dot plug
- Self-closing doors offers added storage protection
- Locking casters for smooth mobility
- Environmentally-friendly natural hydrocarbon refrigerant for reduced global warming potential
- Meets EPA and SNAP requirements
- Convenient security with locks for each section
- Magnetic gaskets to ensure a positive seal
- Pyxis®, Omnicell® and AcuDose RX® compatible



### EXTENDED WARRANTY

- 2-year parts & labor
- 5-year compressor



### DUAL COMPRESSOR

- Independent cooling ensures refrigerated and frozen items are properly stored



### ECO-FRIENDLY

- Uses environmentally friendly hydrocarbons for a greener footprint

# accucold®



SKU	H x W x D	Door Swing	Refrigerator Capacity	Freezer Capacity	Refrigerator Door	Freezer Door	Refrigerator Shelves	Freezer Shelves
ARFG18PV	79.63" x 27.63" x 26.25"	RHD	11.04 cu.ft.	4.5 cu.ft.	Solid with window	Solid	4 wire	2 wire
ARFG15PV	74.63" x 23.5" x 25"	RHD	6.75 cu.ft.	4 cu.ft.	Solid with window	Solid	3 wire	2 wire
ARFG18PVLHD	79.63" x 27.63" x 26.25"	LHD	11.04 cu.ft.	4.5 cu.ft.	Solid with window	Solid	4 wire	2 wire
ARFG15PVLHD	74.63" x 23.5" x 25"	LHD	6.75 cu.ft.	4 cu.ft.	Solid with window	Solid	3 wire	2 wire

## TECHNICAL SPECIFICATIONS (ALL UNITS):

- **Defrost:** Cycle + Manual
- **Thermostat:** Intelligent microprocessor digital
- **Temperature Sensors:** 4 (air probe and buffered probe in each section)
- **Adjustable Operating Control Range:** Refrigerator: from +2°C to +8°C  
Freezer: -30°C to -10°C
- **Temperature Variation:** +/-1°C
- **Shelving:** White plastic coated wire
- **Voltage:** 115V/60Hz
- **Cord:** NEMA 5-15P
- **Lock:** Factory-installed key locks
- **Handle:** Antimicrobial-coated handle to mitigate contamination spread
- **Refrigerant:** Eco-friendly R600a + R290a
- **Light:** LED
- **Access Port:** ½"
- **Door:** Specify RHD/LHD
- **Fan:** Internal
- **Certifications:** ETL, ETL-C
- **Parts & Labor Warranty:** 2 years
- **Compressor Warranty:** 5 years

## INCLUDED

- DO NOT UNPLUG warning sticker
- DO NOT ADJUST TEMPERATURE warning sticker
- DO NOT STORE PERSONAL ITEMS warning sticker

## COMPATIBLE WITH

- Pyxis®
- Omnicell®
- AcuDose RX®

## OPTIONS

- Digital data logger (DDL) installed or as kit with 2-year certificate of conformity
- Wireless Environmental Temperature Monitoring Device (ETM) installed or as kit for cloud-based monitoring and alerts

