


Casco First Responder EMT Stretcher Mattress

The Casco First Responder EMT Stretcher Mattress is made with antimicrobial medical grade foam and a heavy duty tri-laminate vinyl (10oz.) that provides a surface that is hypoallergenic, easy to clean, stain resistant, and durable. The mattress is available in two sizes (73", 78") and includes an Intertek buckle system to provide superior durability and safety for the patient, and ease of use for the EMT crew when moving the patient from the EMT stretcher to another surface. This product is CAL 129 compliant.

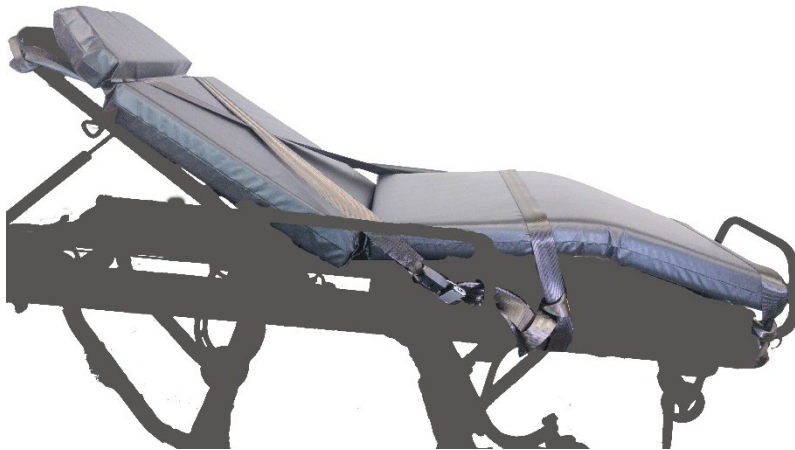

Feature	Advantage	Benefit
Universal Fit	Designed to fit any stretcher frame	Fully adjustable to make, model, and brand of stretcher frame
High density foam core	Comfort	Promotes comfort during transport
Anti-microbial	Protection from all micro-organisms	Helps prevent the spread of germs and further illness
RF Sealed	Protects foam core from fluids	Promotes longer mattress life.
Ventilation Holes	Allows air to escape mattress	Prevents ballooning of mattress
No Bolsters	Easy patient transfer	Lower risk of back injuries to EMTs
Adjustable Buckle Straps	Allows for universal frame fit; can fit any size patient	Allows for standardization of stretcher mattress leading to cost savings

Casco First Responder EMT Stretcher Mattress

	Part #/GTIN	UPC
3x18x73 (mattress only)	H12814/ 00810125103889	810125103889
3x18x73 (mattress + straps)	H12819/ 00810125103919	810125103919
3x18x78 (mattress only)	H12817/ 00810125103902	810125103902
3x18x78 (mattress + straps)	H12816/ 00810125103896	810125103896
Strap Assembly	H12812/ 00810125103872	810125103872



- No bolsters for ease of use



- Adjustable buckle straps allow for product standardization, enabling universal fit and leading to cost savings!



- The Casco First Responder EMT Stretcher mattress is a universal fit mattress that not only fits all stretcher frames, but provides support for patients of all sizes as well.