

CIRCULAR HEAD FLARED EDGE COATED BLOCK

PRODUCT INFORMATION

PRODUCT INFORMATION
CIRCULAR HEAD FLARED EDGE COATED BLOCK
PPFP-PCAK



The Circular Head Flared Edge Coated Block supports and positions the head during medical imaging procedures. It measures 1.5" H, with a 7" top diameter and 9" base diameter. Its circular shape provides a dedicated surface for head positioning.

For added durability, a medical-grade, fluid-resistant coating covers the foam positioner. The coating is easy to clean with medical-grade disinfectant wipes or sprays. It also helps protect the foam interior during repeated clinical use. In addition, the positioner is radiolucent, with low attenuation and minimal shadowing during imaging.

The patented N-VISI-RAY design uses flared edges to significantly reduce artifacts during X-ray imaging. Although the flared sides increase the overall size, they do not affect imaging or positioning.

Finally, the Circular Head Flared Edge Coated Block is MRI Safe and latex-free. This makes it suitable for use across a variety of medical imaging environments.


SPECIFICATIONS

SHIPPING WEIGHT: 16 oz	MATERIAL: Coated Foam
SHIPPING DIMENSIONS: 8 x 8 x 9 in	POSITIONER TYPE: Neck & Head Positioners
PRODUCT DIMENSIONS: 1.5" H, with 7" top diameter and 9" base diameter	FEATURES: Flared Edge, MRI Safe, Latex-Free, Radiolucent, Low Attenuation, Easy to Clean

PRODUCT INFORMATION

CIRCULAR HEAD FLARED EDGE COATED BLOCK

PPFP-PCAK

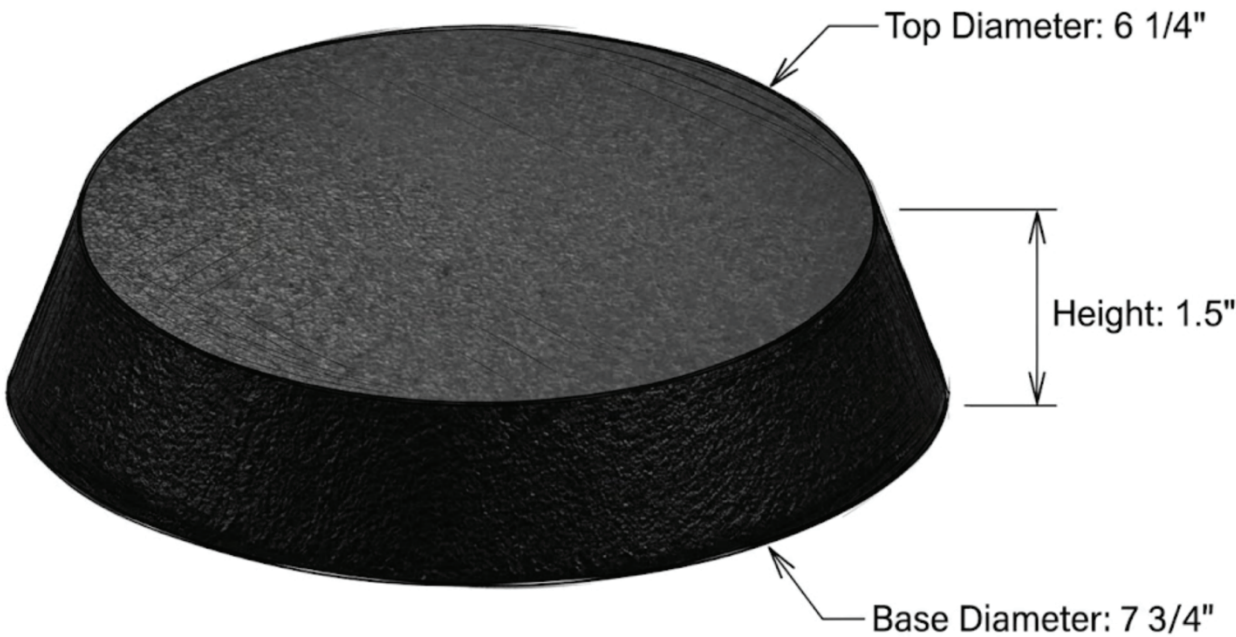

Date: 9/21/2026		Drawing Title: Circular Head Flared Edge Coated Block
Scale: 1:1	Unit: Inches	Drawing No: PPFP-PCAK
		Notes: All dimensions are ±0.1"