


ANGULAR FLARED EDGE COATED BOLSTER

PRODUCT INFORMATION

PRODUCT INFORMATION
ANGULAR FLARED EDGE COATED BOLSTER
PPFP-PCAY

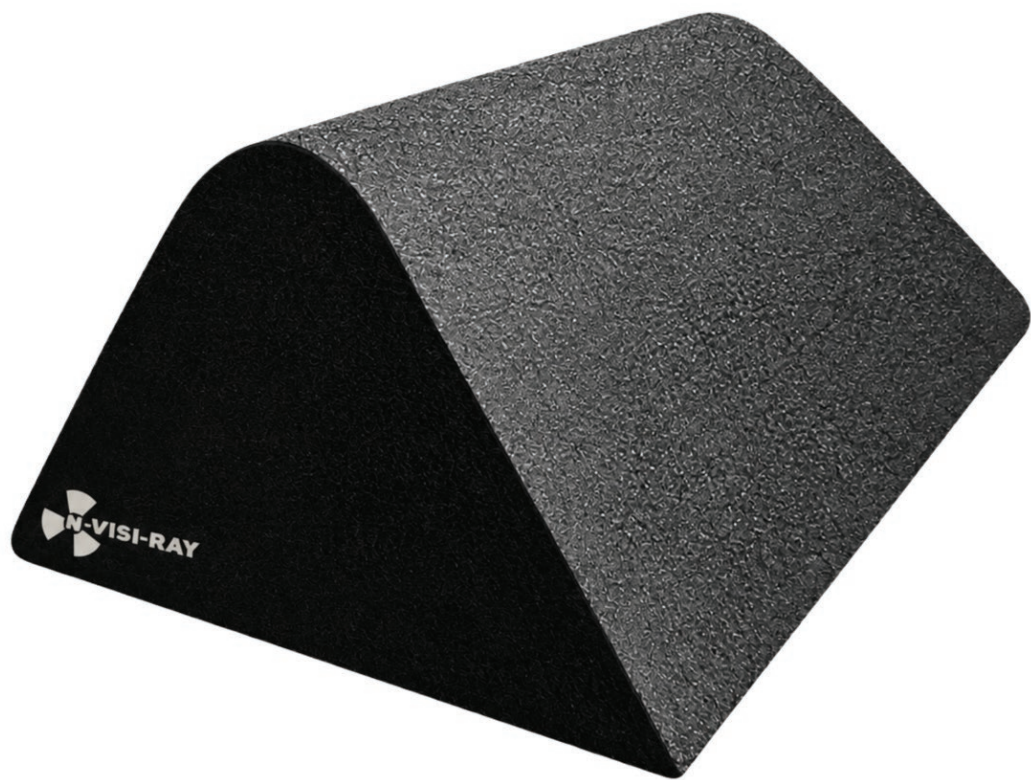

The Angular Flared Edge Coated Bolster supports and positions the lumbar and spinal areas during medical imaging procedures. Measuring 22.5" L x 14" W x 7.25" H, its angled design provides support for applications such as lumbar imaging and spinal injections.

For added durability, a medical-grade, fluid-resistant coating covers the foam positioner. The coating is easy to clean with medical-grade disinfectant wipes or sprays. It also helps protect the foam interior during repeated clinical use. In addition, the positioner is radiolucent, with low attenuation and minimal shadowing during imaging.

The patented N-VISI-RAY design uses flared edges to significantly reduce artifacts during X-ray imaging. Although the flared sides increase the overall size, they do not affect imaging or positioning.

Finally, the Angular Flared Edge Coated Bolster is MRI Safe and latex-free. This makes it suitable for use across a variety of medical imaging environments.

SPECIFICATIONS


SHIPPING WEIGHT: 64 oz
SHIPPING DIMENSIONS: 26 x 14 x 11 in
PRODUCT DIMENSIONS: 22.5" L x 14" W x 7.25" H

MATERIAL: Coated Foam
POSITIONER TYPE: Torso Positioners
FEATURES: Flared Edge, MRI Safe, Latex-Free, Radiolucent, Low Attenuation, Easy to Clean

PRODUCT INFORMATION

ANGULAR FLARED EDGE COATED BOLSTER

PPFP-PCAY

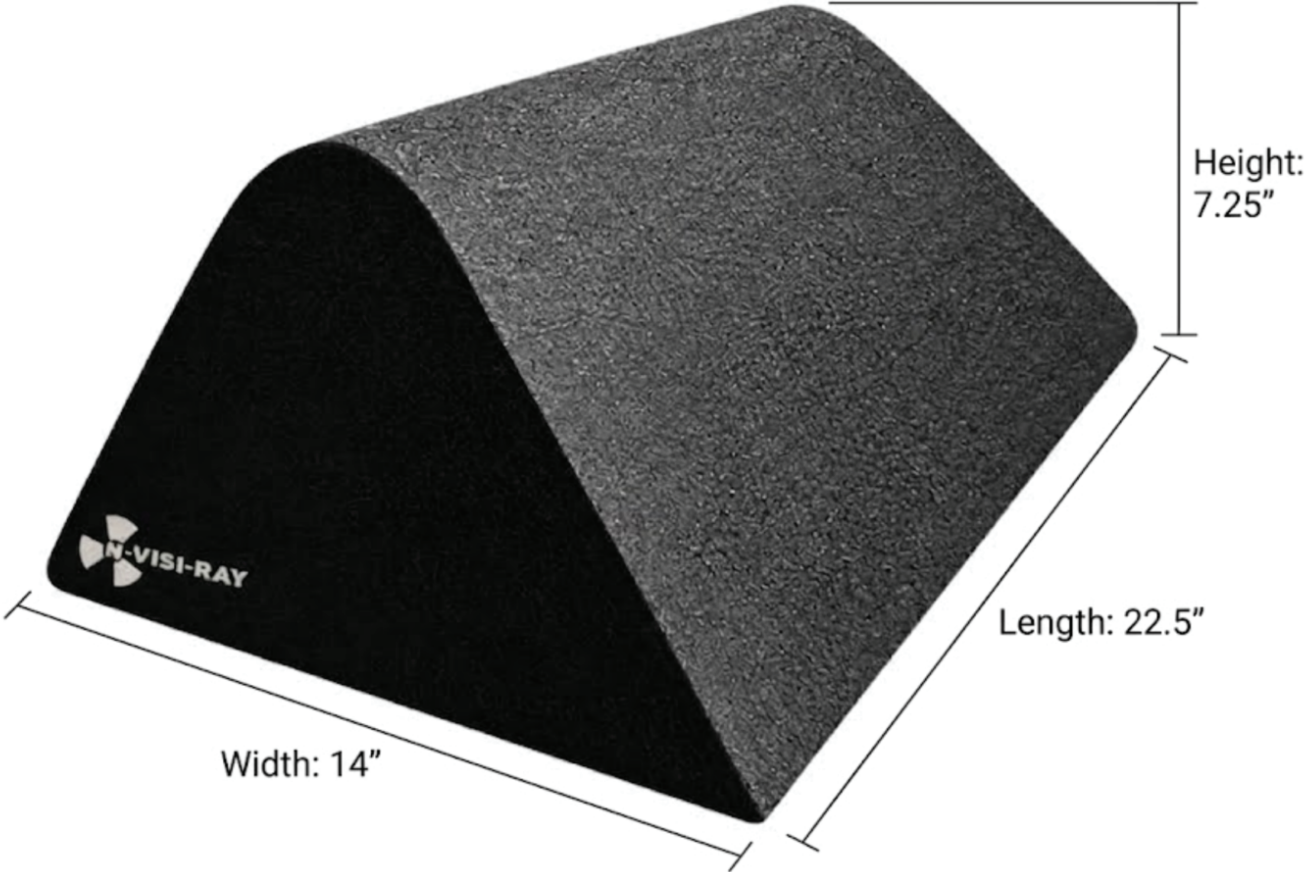

Date: 9/21/2026		Drawing Title: Angular Flared Edge Coated Bolster
Scale: 1:1	Unit: Inches	Drawing No: PPFP-PCAY
		Notes: All dimensions are ± 0.1 "