Quick Reference Guide

Team3 USA

Fetal Monitor

Switching the unit ON

1 Touch and hold the On/Off button for 2 seconds. (On Team 3 Front Panel).

2 Check the printer: Is there sufficient paper for the monitoring session?

3 Connect transducers. Display shows transducer connection

Control Bar



Printer On/ Off	Paper Feed	Numbers or Graph	Volume Down	Volume Up	Notes *	Settings	Menu *

(Depending on model / operating mode)*

1 Touch any of the keys above to operate that function.

Entering Patient Details



1 Touch and hold the box in the top left corner of the screen and a Patient Details dialogue appears.

2 Enter Patient Name using the popup keyboard.

3 Touch the Patient ID field and type data as required.

4 Touch the Gestational Age field to bring up another dialogue. Use the up and down arrows to set the correct gestational age.

5 Touch ☑ to set Gestational Age.

6 Touch ☑ to set Patient Data.

7 Ensure patient data is entered before you print. It cannot be changed afterwards.

Monitoring FHR

1 Determine fetal lie and position.

2 Apply Gel to the green transducer face.


Position on abdomen so as to get optimal heart sound. Secure with belt.

3 Make sure the FHR is clear, and different from the maternal pulse taken at the mother's wrist.

FHR signal quality is shown by the flashing heart symbol.

4 Adjust the volume for the audio channel.

Monitoring Uterine Activity

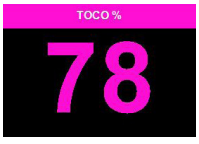

1 Position the pink transducer on the midline over the fundus of the uterus. Secure with belt.

2 Zero the Toco. Make sure the uterus is relaxed, and then touch and hold the screen in the pink TOCO Clinical Region.

HUNTLEIGH

A MEMBER OF THE ARJO FAMILY

Monitoring FHR for twins

1 Connect both green transducers to Team 3. Display shows transducer connection.


2 Apply Gel to the transducer faces, and position them so as to get optimal heart sounds.

3 Ensure you are not getting two signals from the same twin. If you are, Team 3 will warn that similar FHRs have been detected.

4 Secure the transducers with belts.

5 To change the Audio channel, touch and hold the screen in the FHR Display Region of relevant trace.

Loading paper

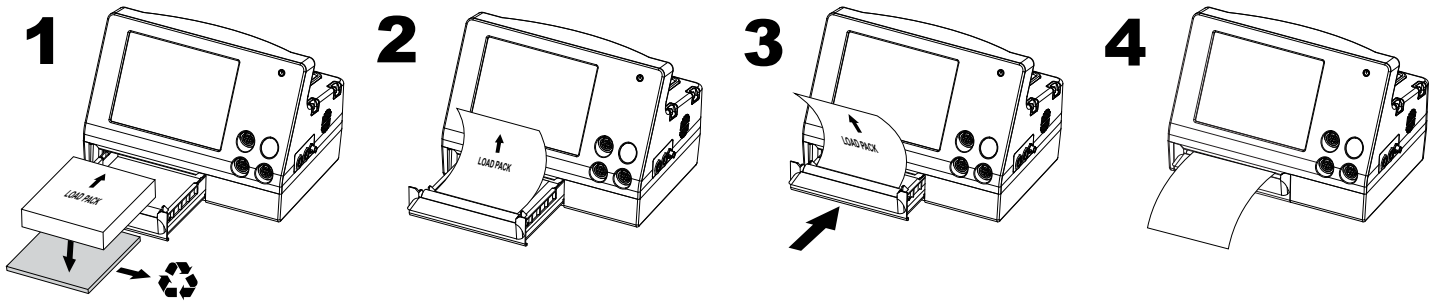

Printing the trace

1 To start or stop printing *, touch Printer On/Off key  on the control bar.

2 To advance the paper, touch Paper Feed .

** Refer to the Settings Menu in the instructions for Use (Settings > System> Recording Limits) for details on setting a automatic printing stop time.*

Store Trace (PDF) to USB Drive (8GB max) *


1 Select View Menu  .

2 Review Trace  .


3 Enter PIN code to enable access to stored traces (PIN code : 98765).

4 Select trace from list by using the arrow to scroll through the traces  .

5 Insert USB into USB drive in the rear panel of the fetal monitor.

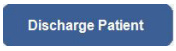6 Select View Menu and then select Discharge Patient  .

7 Remove USB from the fetal monitor to enable trace to be shared on a PC.

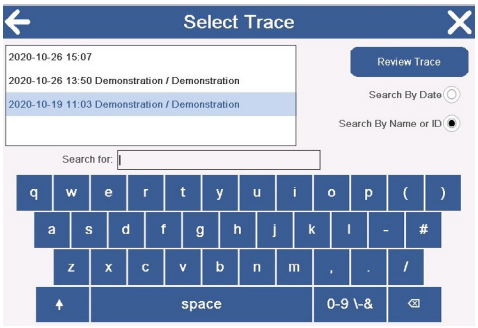

** Refer to Sonicaid Reporter Quick Reference Guide for details on uploading stored traces to the Sonicaid Reporter software package.*

Switching off

1 Touch and hold the On/Off button, (On Team 3 Front Panel)  , until ‘Power’ screen is displayed. Touch ‘Power Off’